

2019 PRODUCT CATALOG



INTRODUCTION

About Garreco	1
Sales Policy	2
Custom Products	3
Product Samples	3
Online Resources	3

GYPSUM

Die Stones	4
Apex™, Type IV-R Resin	4
Cad-Scan™, Type IV-R Resin	5
Excalibur™, Type IV High Strength	5
Express Set™, Type IV Specialty	6
NaturalRock EH™, Type IV-R Resin	6
Satinstone™, Type IV	6
Veri-Die™, Type IV	7
Veri-Die™ 20, Type V High Expansion	7
Base Stones	8
Tecstone™, Type III	8
Tecstone™ FL - Flowable, Type III	8
Fastone™, Type III	9
Labstone™, Type III	9
Flask It™, Type III	10
Flasking Stone™, Type III	10
Mounting Stone™, Type III	10
Orthodontic Gypsum	11
Orthoplaster™, Type III	11
Orthostone™, Type III	11
Specialty Gypsum	11
X-Roc™, Expansion Stone	11
Plaster	12
Plaster, Type II	12
Articulating Plaster, Type II	12
Operatory Gypsum	13
Pristine® II, Type II	13
Pristine® III, Type III	13
Pristine® V, Type IV	13

ABRASIVES

Resilience™ Acrylic Polishing Media	14
Aluminum Oxide	15
Pumice	16
Glass Beads	17
Walnut Shells	17

ACRYLICS

Dynasty™ Hi-Impact Acrylic	18
Hi-Impact Denture Base	18
Hi-Impact Pour Type Denture Base	19
Liquid Monomer	19
Axcent™ Acrylic	20
Premium Denture Base	20
Standard Denture Base	20
15 Minute Denture Base	21
Repair Denture Base	21
Pour Type Denture Base	22
Orthodontic Acrylic Powder	22
Tray Material	23
Softliner Reline Material Kit	23
Liquid Monomer	24

CROWN & BRIDGE

Investment Products	25
Ultravest™, Investment Material	25
Perfect Press™, Investment Material	25
Investment Liquid	26
Flexi-Ring™ Ringless C&B System	26
C&B Supplies	27
Die Spacer	27
Debubblizer	27
Plaster Separator	27
Vest-It™ Refractory Material	28
Honeycomb Trays	28
Flex Tray	28
Disposable Plungers	28
Crucibles	29
Natural Agate	29
Phyx-It Model and Die Adhesive	29

LAB SUPPLIES

General Use	30
Separating Film	30
PERform™ Products	31
Tooth Sieve	31
Tooth Cutting Bur	31
Filter Papers	31

PARTIAL DENTURES

Vitex™ Products	32
Investment Material	32
3D Vest Investment Material	33
Refractory Model Hardener	33
Duplicating Material	33
Alloy	34
Superforms	34

ROTARY

Perlada (PD) Polishers	35
NaturalGlaze Diamond Polishing System	36
Interproximal Rotary	36

SILICONE

Silaform™ Silicone Putty	37
Accusil™ Silicone Duplicating Material	38
Rapidsil™ Silicone Duplicating Material	38
Silsaver™	39
Duplicating Cross	39
Sprue Formers	39
Duplicating Ring	40
25KG Tote Stand	40
Silicone Dispenser	40
Silicone Duplication Technique	40

RESOURCES

Find a Dealer	41
Contact Us	41
Referral Program	42



ABOUT GARRECO

Garreco Dental has been providing quality dental lab products since 1988. As one of the leading manufacturers of dental gypsum and abrasive products, we take great pride in the quality and consistency of all of our products. Our manufacturing process is ISO 13485 certified and uses only the best raw materials which we carefully formulate for exacting product standards and superior performance. In addition, 100% of our gypsum products are tested prior to being shipped to ensure quality and consistency.



In addition to our core dental business, Garreco also produces products for the fields of jewelry casting, microdermal abrasion, forensics, taxidermy, special effects, and hobby. Please contact our sales team for more information regarding these products.

Our dedicated team of professionals is committed to providing superior customer service and on-time delivery for all our products. We encourage you to reach out and allow us to impress you with our products, honest approach to business, and the Garreco Advantage.



THE GARRECO ADVANTAGE

- ★ Quality Dental Lab Products since 1988
- ★ ISO Certified Manufacturer
- ★ Highest Quality Product for a Great Price
- ★ Consistency in Every Box
- ★ Right Product Delivered on Time
- ★ Custom Blending

SALES POLICY



TERMS

Our normal terms are Net-30 days. However, sales will be handled on a COD/Prepay basis until an accepted credit application is on file. See Pricing for additional information.

MINIMUM ORDER

Our minimum order amount is \$50.00 for combined products per invoice, excluding freight charges. A \$5.00 handling fee will be assessed on all orders below stated minimum.

PRICING

We make every attempt to reduce costs and keep our prices as low as possible. However, we find that incremental increases are necessary. Several factors are taken into account such as: raw material and packaging cost, transportation cost, volume per item, volume per order, and general cost of doing business.

FREIGHT CHARGES

Because freight costs are often a large percentage of your total order, we endeavor to obtain the lowest shipping cost possible. Orders will be shipped, freight prepaid and all charges will be added to invoice, unless prior arrangements have been made.

PALLET FEE

A pallet fee will be charged on all palletized orders.

DAMAGE CLAIMS

Any damaged goods should be immediately noted on delivery receipt when freight is received (for FedEx deliveries, please notify driver). If damage is not noted, Garreco will not be responsible for corrective action. Damage claims must be submitted to Garreco within 7 days of delivery. Furthermore, Garreco will not be responsible for corrective action for carriers other than those preferred by Garreco.

PRODUCT RETURNS

Prior to returning any merchandise to Garreco, an RA Number must be obtained from our office and affixed or noted on the outer case for all returned goods. Please note, Garreco charges 15-20% restocking fees for all returned items depending on product labeling. Any returned items received by Garreco in unsatisfactory condition will be rejected.

CUSTOM PRODUCTS

See what Garreco Dental can do for you! As a private label manufacturer for some of the largest dental companies in the world, Garreco is well-versed in offering a premium product customized to your specific needs. We offer an in-house design team to design your private label artwork.



CUSTOM COLOR

Our Custom Color Program offers a choice from a wide variety of pastel shades for their gypsum products, including matching company colors.



CUSTOM SET

Custom Set Products are available for almost all Garreco gypsums. This service offers a specified set time from 2 to 40 minutes. Whether you need to finish one case quickly or a hundred cases at the same time, we can customize a product for you.



CUSTOM PACKAGING

Custom Packaging services are available for all Garreco products.

PRODUCT SAMPLES



Samples for most of our products are available upon request. Please note that some products are not available to sample. Contact us for more information.

To make your sample request or get more information regarding samples, please contact the Garreco Sales Team at 800-334-1443 ext. 806 or email us at sales@garreco.com.

ONLINE RESOURCES

GARRECO.COM

Visit us at Garreco.com to view our Featured Products and get support. We have Safety Data Sheets available to download for each product.

PROMOTIONS & SPECIAL OFFERS

Visit us at Garreco.com/promotions to view our current specials on the dental products you love.

GARRECO DENTAL BLOG

Join us at Garreco.com/blog to read about what's happening at Garreco and get involved in the discussion!



LIKE, FOLLOW, & SHARE
GARRECO DENTAL ON:





DIE STONES

APEX™
Type IV-R Resin Reinforced Die Stone

Apex is a resin reinforced natural gypsum die material, available in two expansion variations (RZ 10, RZ 15) making it the perfect choice to complement any technique. Perfect for all crown and bridge applications as well as implant or precision attachment cases. Resin additives eliminate brittleness and add superior external hardness. The exclusive Apex flow agent allows for excellent stackability, while promoting free material flow during vibration.

INDICATION FOR USE: Crown and bridge; porcelain; implants; cast partials; high-precision cases

SETTING TIME: 10-13 min
WATER/POWDER RATIO: 19ml:100g
COMPRESSIVE STRENGTH: 18,500 psi
EXPANSION: 0.10% (RZ 10), 0.15% (RZ 15)



RZ 10 COLORS AND SIZES AVAILABLE:

	BLUE	GREEN	PINK	IVORY*	WHITE	TAN
PACKAGING						
120X70G FOIL POUCH	1301580	1302580	1305080	1303580	1303080	1306580
25LB CARTON	1301525	1302525	1305025	1303525	1303025	1306525
50LB CARTON	1301550	1302550	1305050	1303550	1303050	1306550

RZ 15 COLORS AND SIZES AVAILABLE:

	BLUE	GREEN	GREY	IVORY*	WHITE	TAN
PACKAGING						
120X70G FOIL POUCH	1321580	1322580	1324580	1323580	1323080	1326580
25LB CARTON	1321525	1322525	1324525	1323525	1323025	1326525
50LB CARTON	1321550	1322550	1324550	1323550	1323050	1326550

WARNING: Cancer - www.P65Warnings.ca.gov

FOR SALES CALL 800-334-1443 EXT. 806
 OR EMAIL US AT SALES@GARRECO.COM

DIE STONES

CAD-SCAN™

Type IV-R Resin Die Stone

CAD-scan is a unique gypsum material that contains a special scanning additive and resin reinforcement. CAD-scan reflects light nicely for essential readability, allowing for a shorter working time, and is manufactured from an ultra-fine gypsum base, which prevents nicks and abrasions caused by brittleness. CAD-scan has a low expansion and high compressive strength, making it the material of choice for all CAD/CAM restorations. It is perfect for production environments, as well as those single cases requiring a rapid set.

INDICATION FOR USE: Crown and bridge; porcelain; implants; cast partials; high-precision cases

SETTING TIME: 11-13 min

WATER/POWDER RATIO: 20ml:100g

COMPRESSIVE STRENGTH: 18,500 psi

EXPANSION: 0.06%-0.09%

COLORS AND SIZES AVAILABLE:



PACKAGING	GOLDEN	GREEN	BUFF
10LB CANISTER	1466710	1462510	1461010
25LB PAIL	1466725	1462525	1461025
50LB CARTON	1466750	1462550	1461050

⚠ WARNING: Cancer - www.P65Warnings.ca.gov

EXCALIBUR™

Type IV High Strength Die Stone

With its smooth mixing characteristics and ample setting time, Excalibur produces a “textbook” model every time. Die trimming becomes a worry-free procedure due to the ultra-fine gypsum base, which prevents nicks and abrasions caused by brittleness. Recommended for use as crown and bridge or partial denture models. Suitable for use with all impression materials. Gypsum hardeners are not recommended.

INDICATION FOR USE: Crown and bridge; porcelain; implants; cast partials; high-precision cases

SETTING TIME: 12-15 min

WATER/POWDER RATIO: 22ml:100g

COMPRESSIVE STRENGTH: 18,000 psi

EXPANSION: 0.09%

COLORS AND SIZES AVAILABLE:



PACKAGING	BLUE	GREEN	VIOLET	GREY	WHITE	TAN
120X70G FOIL POUCH	1351580	1352580	1352080	1354580	1353080	1356580
25LB CARTON	1351525	1352525	1352025	1354525	1353025	1356525
50LB CARTON	1351550	1352550	1352050	1354550	1353050	1356550

DIE STONES

EXPRESS SET™

Type IV Specialty Stone

ExpressSet is a Type IV specialty stone with an extremely fast set time, smooth mixing consistency, and uniform color matrix. Models can be released from the impression in 5 minutes, producing a precision cast expressly suited for mouth guards, custom trays, study casts, models for opto-electronic/laser scanning with powder, temporaries and denture repairs.



INDICATION FOR USE: Crown and bridge; porcelain; implants; cast partials; high-precision cases

COLORS AND SIZES AVAILABLE:

SETTING TIME: 60-90 seconds
WATER/POWDER RATIO: 23ml:100g
COMPRESSIVE STRENGTH: 14,000 psi
EXPANSION: 0.15%

PACKAGING	BUFF	PINK
80X140G FOIL POUCH	1401580	1402580
25LB CARTON	1401525	1402525
50lb CARTON	1401550	1402550

NATURALROCK™ EH

Type IV-R Resin Reinforced Die Stone

NaturalRock EH is a traditional, creamy, smooth, resin reinforced die stone which exhibits high compressive strength and very low expansion.



INDICATION FOR USE: Crown and bridge; porcelain; implants; cast partials; high-precision cases

COLOR AND SIZE AVAILABLE:

SETTING TIME: 10-13 min
WATER/POWDER RATIO: 20ml:100g
COMPRESSIVE STRENGTH: 18,000 psi
EXPANSION: 0.08%

PACKAGING	LT. TAN
120X70G FOIL POUCH	1446780
25LB CARTON	1446725
50LB CARTON	1446750

SATINSTONE™

Type IV Die Stone

Comparable to Whip Mix's Silky-Rock, SatinStone™ is distinguished by exceptional toughness and an extraordinary "satin-smooth" finish that works well with all types of impression materials. Introducing only the finest raw materials allows this product to have a creamy and effortless mix each and every time. SatinStone™ also has a low expansion which makes it one of the most the most unique die stones on the market.



INDICATION FOR USE: Crown and bridge; porcelain; implants; cast partials; high-precision cases

COLORS AND SIZES AVAILABLE:

SETTING TIME: ~10 min
WATER/POWDER RATIO: 23ml:100g
COMPRESSIVE STRENGTH: 15,000 psi
EXPANSION: 0.09%-0.11%

PACKAGING	ROSE	WHITE
25LB CARTON	1175025	1173025
50LB CARTON	1175050	1173050

DIE STONES

VERI-DIE™
Type IV Die Stone

Produced to handle the many different laboratory techniques that are being used in today’s dental labs. With its range in application and versatility, Veri-Die is perfect for almost every requirement. Its physical properties allow Veri-Die to conform to the requirements for each use. Due to its range in particle size, Veri-Die virtually eliminates normal chipping and abrasion. Veri-Die has become the worry free material of choice by providing the technician ample working time and a moderate expansion rate.


INDICATION FOR USE: Crown and bridge;
 porcelain; implants; cast partials;
 high-precision cases

SETTING TIME: 10-13 min
WATER/POWDER RATIO: 21ml:100g
COMPRESSIVE STRENGTH: 16,500 psi
EXPANSION: 0.15%

COLORS AND SIZES AVAILABLE:



PACKAGING	BLUE	GREEN	PINK	IVORY*
120X70G FOIL POUCH	1401580	1402580	1405080	1403580
25LB CARTON	1401525	1402525	1405025	1403525
50LB CARTON	1401550	1402550	1405050	1403550

*  **WARNING:** Cancer - www.P65Warnings.ca.gov

VERI-DIE™ 20
Type V High Expansion Die Stone

Produced to handle the many different laboratory techniques that are being used in today’s dental labs. With its range in application and versatility, Veri-Die is perfect for almost every requirement. Its physical properties allow Veri-Die to conform to the requirements for each use. Due to its range in particle size, Veri-Die virtually eliminates normal chipping and abrasion. Veri-Die has become the worry free material of choice by providing the technician ample working time and a moderate expansion rate.

INDICATION FOR USE: Crown and bridge;
 porcelain; implants; cast partials;
 high-precision cases

SETTING TIME: 10-13 min
WATER/POWDER RATIO: 21ml:100g
COMPRESSIVE STRENGTH: 18,000 psi
EXPANSION: 0.20%

COLOR AND SIZE AVAILABLE:



PACKAGING	GREEN
120X70G FOIL POUCH	1432580
25LB CARTON	1432525
50LB CARTON	1432550

BASE STONES

TECSTONE™

Type III Base Stone

Widely recognized as one of the best stones on the market and designed for ultimate versatility, Tecstone is recommended for use in models and flasking, in the production of full or partial dentures, and in basing crown and bridge models. A creamy mix, moderate compressive strength, and superior outward hardness are all characteristics of this unique formula.

INDICATION FOR USE: Full or partial denture models; flasking; basing crown and bridge models

SETTING TIME: 13-15 min

WATER/POWDER RATIO: 30ml:100g

COMPRESSIVE STRENGTH: 9,000 psi

EXPANSION: 0.12%

COLORS AND SIZES AVAILABLE:

PACKAGING	BUFF	WHITE	BLUE
25LB CARTON	1011025	1013025	1011525
50LB CARTON	1011050	1013050	1011550



TECSTONE™ FL - FLOWABLE

Type III Base Stone

This unique new product from Garreco is specifically designed for basing crown & bridge models. Tecstone FL - Flowable pours into base formers without vibration and yields a bubble-free base. Tecstone FL - Flowable far exceeds the strength of conventional basestones while offering a dense, smooth result. Accuracy is also addressed by a low expansion.

INDICATION FOR USE: Basing crown and bridge models

SETTING TIME: 9-11 min


WATER/POWDER RATIO: 23ml:100g


COMPRESSIVE STRENGTH: 12,500 psi

EXPANSION: 0.08%

COLORS AND SIZES AVAILABLE:

PACKAGING	BUFF	WHITE*	BLUE
25LB CARTON	1031025	1033025	1031525
50LB CARTON	1031050	1033050	1031550

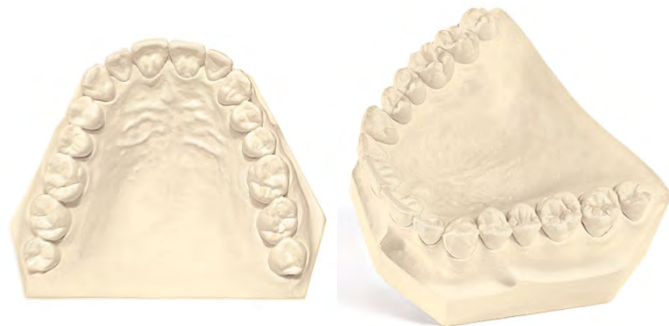

*  **WARNING:** Cancer - www.P65Warnings.ca.gov

BASE STONES

FASTONE™

Type III Base Stone

This dental stone offers the same unique characteristics as our Tecstone product, with a faster setting time. Used for full or partial denture models, denture flasking procedures, and recommended for use in basing crown and bridge models.

INDICATION FOR USE: Full or partial denture models; flasking; basing crown and bridge models

SETTING TIME: 10-12 min

WATER/POWDER RATIO: 30ml:100g

COMPRESSIVE STRENGTH: 9,000 psi

EXPANSION: 0.12%

COLORS AND SIZES AVAILABLE:

PACKAGING	BUFF	GOLDEN	BLUE	WHITE
25LB CARTON	1101025	1106725	1101525	1103025
50LB CARTON	1101050	1106750	1101550	1103050



LABSTONE™

Type III Base Stone

The “work horse” of the gypsum family. Designed for ultimate versatility. Labstone is recommended for use in models and flasking, in the production of full or partial dentures, and in basing crown and bridge models. Creamy mix, moderate compressive strength, and superior outward hardness, are all characteristics of this unique formula.

INDICATION FOR USE: Full or partial denture models, denture flasking procedures, basing crown and bridge models

SETTING TIME: 11-14 min

WATER/POWDER RATIO: 30ml:100g

COMPRESSIVE STRENGTH: 8,500 psi

EXPANSION: 0.12%

COLORS AND SIZES AVAILABLE:

PACKAGING	BUFF	WHITE	BLUE
25LB CARTON	1051025	1053025	1051525
50LB CARTON	1051050	1053050	1051525



PLASTER

PLASTER

Type II Gypsum

Dental Plasters offered by Garreco are formulated especially for use in today’s complex field of dentistry. Designed for articulating and flasking procedures, Garreco plasters mix smoothly and exhibit excellent working characteristics. Available in three formulations with varying set times.

INDICATION FOR USE: Articulating and flasking procedures

SETTING TIME: 17-19 min (Slow); 13-15 min (Medium); 9-11 min (Fast)

WATER/POWDER RATIO: 45ml:100g

COMPRESSIVE STRENGTH: 3,500 psi

EXPANSION: 0.19%

COLORS AND SIZES AVAILABLE:



PACKAGING	SLOW - WHITE	REGULAR - WHITE	FAST - WHITE
25LB CARTON	1550025	1450025	1500025
50LB CARTON	1550050	1450050	1500050

ARTICULATING PLASTER

Type II Gypsum

Articulating Plaster exhibits similar characteristics to our Dental Plasters, but with a greatly reduced set time. Recommended for articulation use.

INDICATION FOR USE: Articulation

SETTING TIME: 3-5 min

WATER/POWDER RATIO: 45ml:100g

COMPRESSIVE STRENGTH: 3,500 psi

EXPANSION: 0.19%

COLOR AND SIZE AVAILABLE:



PACKAGING	WHITE
25LB CARTON	1560025
50LB CARTON	1560050

ORTHODONTIC GYPSUM

ORTHOPLASTER

Type III Base Stone

High strength plaster designed for orthodontic applications. Orthoplaster combines smooth mixing, low expansion and ample setting time. Contains special whitening additives for a brilliant white color.



INDICATION FOR USE: Orthodontic applications

SETTING TIME: 10-13 min (Regular); 6-9 (Fast)

WATER/POWDER RATIO: 35ml:100g

COMPRESSIVE STRENGTH: 6,500 psi

EXPANSION: 0.12%

COLOR AND SIZES AVAILABLE:

PACKAGING	REGULAR	FAST
25LB CARTON	1600025	1620025
50LB CARTON	1600050	1620050



WARNING: Cancer - www.P65Warnings.ca.gov

ORTHOSTONE

Type III Base Stone

OrthoStone exhibits the strength and accuracy of a quality base stone, plus special additives resulting in a brilliant white color. The superior color, working properties, and hardness make this the product of choice for those specializing in orthodontics.



INDICATION FOR USE: Orthodontic applications

SETTING TIME: 13-15 min

WATER/POWDER RATIO: 30ml:100g

COMPRESSIVE STRENGTH: 9,000 psi

EXPANSION: 0.12%

COLORS AND SIZES AVAILABLE:

PACKAGING	WHITE
25LB CARTON	1253025
50LB CARTON	1253050



WARNING: Cancer - www.P65Warnings.ca.gov

SPECIALTY GYPSUM

X-ROC™

Expansion Stone

X-Roc was developed for pouring the duplicate model in processing thermoplastic appliances such as: Flexite®, Valplast®, Flex-Star V, TCS®, FRS®, DuraFlex®, & DurAcetal®. This stone compensates for the shrinkage of these materials, thereby eliminating tight fits on removable prosthetics processed with conventional stones.



INDICATION FOR USE: Specialty duplication material

SETTING TIME: 22-25 min

WATER/POWDER RATIO: 29ml:100g

COMPRESSIVE STRENGTH: 13,000 psi

EXPANSION: 1.10%

COLORS AND SIZES AVAILABLE:

PACKAGING	WHITE
25LB CARTON	1483025
50LB CARTON	1483050

PLASTER

PLASTER

Type II Gypsum

Dental Plasters offered by Garreco are formulated especially for use in today's complex field of dentistry. Designed for articulating and flasking procedures, Garreco plasters mix smoothly and exhibit excellent working characteristics. Available in three formulations with varying set times.

INDICATION FOR USE: Articulating and flasking procedures

SETTING TIME: 17-19 min (Slow);
13-15 min (Medium); 9-11 min (Fast)

WATER/POWDER RATIO: 45ml:100g

COMPRESSIVE STRENGTH: 3,500 psi

EXPANSION: 0.19%

COLORS AND SIZES AVAILABLE:



PACKAGING	SLOW - WHITE	REGULAR - WHITE	FAST - WHITE
25LB CARTON	1550025	1450025	1500025
50LB CARTON	1550050	1450050	1500050

ARTICULATING PLASTER

Type II Gypsum

Articulating Plaster exhibits similar characteristics to our Dental Plasters, but with a greatly reduced set time. Recommended for articulation use.

INDICATION FOR USE: Articulation

SETTING TIME: 3-5 min

WATER/POWDER RATIO: 45ml:100g

COMPRESSIVE STRENGTH: 3,500 psi

EXPANSION: 0.19%

COLOR AND SIZE AVAILABLE:



PACKAGING	WHITE
25LB CARTON	1560025
50LB CARTON	1560050

OPERATORY GYPSUM

PRISTINE® II Type II Gypsum

A fast setting Type II model and articulating plaster exhibiting excellent working characteristics. Has a fresh, clean scent to mask odors picked up by impression materials which can be transferred to working models. Comes in stackable pails with flip top re-closable lid and a convenient scoop.



INDICATION FOR USE: Gypsum designed for dental office applications

SETTING TIME: 6-10 min
WATER/POWDER RATIO: 45ml:100g
COMPRESSIVE STRENGTH: 3,500 psi
EXPANSION: 0.19%

COLOR AND SIZE AVAILABLE:

PACKAGING	WHITE
20LB PAIL	1500020

PRISTINE® III Type III Base Stone

With a broad range of applications this Type III stone may be used for full and partial denture models, denture flasking procedures, and for basing crown and bridge models. Has a fresh, clean scent to mask odors picked up by impression materials which can be transferred to working models. Comes in stackable pails with flip top re-closable lid and a convenient scoop.



INDICATION FOR USE: Gypsum designed for dental office applications

SETTING TIME: 8-10 min
WATER/POWDER RATIO: 30ml:100g
COMPRESSIVE STRENGTH: 9,000 psi
EXPANSION: 0.12%

COLORS AND SIZES AVAILABLE:

PACKAGING	GOLDEN	WHITE
25LB PAIL	1026725	1023025

PRISTINE® V Type IV Die Stone

Versatility for all your fixed and removable models. This Type IV die stone resists normal chipping and abrasion. Has a fresh, clean scent to mask odors picked up by impression materials which can be transferred to working models. Comes in stackable pails with flip top re-closable lid and a convenient scoop.



INDICATION FOR USE: Gypsum designed for dental office applications

SETTING TIME: 10-13 min
WATER/POWDER RATIO: 21ml:100g
COMPRESSIVE STRENGTH: 15,000 psi
EXPANSION: 0.15%

COLORS AND SIZES AVAILABLE:

PACKAGING	GREEN	PINK
25LB PAIL	1382525	1385025



ABRASIVES

RESILIENCE™
Polishing/Finishing Media

Resilience™ is an all-purpose acrylic polishing media designed for faster and smoother polishing of acrylic appliances. Specially formulated with unique abrasives to cut faster, easier, and smoother than pumice-type products with no scratching or grooving. Resilience saves time and does not require a change in your current technique. Kind to your hands and does not give off a bad odor commonly associated with pumice.

INDICATION FOR USE: Polishing acrylic appliances

COLORS AND SIZES AVAILABLE:

ITEM #	SIZE	PACKAGING
2019003	2.5LB	PAIL
2019010	10LB	CARTON
2019020	20LB	CARTON

BENEFITS OF USING RESILIENCE:

- It's easy to use
- Replacement for pumice media
- One grit size does it all
- No odor or bacteria
- No scratching or grooving
- Saves time
- High shine polishing



 **WARNING:** Cancer - www.P65Warnings.ca.gov








ABRASIVES

ALUMINUM OXIDE

Garreco Aluminum Oxide blast abrasives are among the finest in the world and used for divesting and cleaning of appliances. Our 99.6% pure (white) and 96% pure (brown) electro-fused Aluminum Oxide abrasives are machine-sized and computer-scanned to provide uniform particle size. Available in popular sizes and packed with desiccant to prevent moisture absorption.

INDICATION FOR USE: Divesting/cleaning of appliances

COLORS AND SIZES AVAILABLE:

ITEM #	COLOR	GRIT	SIZE	PACKAGING
2101025	WHITE	25 MICRON	25LB	CARTON
2151025	WHITE	25 MICRON	25LB	PAIL
2101050 	WHITE	25 MICRON	50LB	CARTON
2151050	WHITE	25 MICRON	50LB	PAIL
2172002S	WHITE	50 MICRON	2LB	JAR
2172002	WHITE	50 MICRON	6X2LB	JAR
2102025	WHITE	50 MICRON	25LB	CARTON
2152025	WHITE	50 MICRON	25LB	PAIL
2102050 	WHITE	50 MICRON	50LB	CARTON
2152050	WHITE	50 MICRON	50LB	PAIL
2103025	WHITE	60 MICRON	25LB	CARTON
2153025	WHITE	60 MICRON	25LB	PAIL
2103050 	WHITE	60 MICRON	50LB	CARTON
2153050	WHITE	60 MICRON	50LB	PAIL
2104025	WHITE	100 MICRON	25LB	CARTON
2154025	WHITE	100 MICRON	25LB	PAIL
2104050 	WHITE	100 MICRON	50LB	CARTON
2154050	WHITE	100 MICRON	50LB	PAIL
2104525	WHITE	110 MICRON	25LB	CARTON
2154525	WHITE	110 MICRON	25LB	PAIL
2104550 	WHITE	110 MICRON	50LB	CARTON
2154550	WHITE	110 MICRON	50LB	PAIL
2106025	WHITE	250 MICRON	25LB	CARTON
2156025	WHITE	250 MICRON	25LB	PAIL
2106050 	WHITE	250 MICRON	50LB	CARTON
2156050	WHITE	250 MICRON	50LB	PAIL
2202025	BROWN	50 MICRON	25LB	CARTON
2252025	BROWN	50 MICRON	25LB	PAIL
2202050 	BROWN	50 MICRON	50LB	CARTON
2252050	BROWN	50 MICRON	50LB	PAIL

 Denotes an item that cannot be shipped by FedEx Ground; must be shipped FedEx Freight only



ABRASIVES

PUMICE

A consistently fast, smooth cutting, long lasting abrasive. Suitable for polishing full or partial dentures. Available in a variety of particle sizes.

INDICATION FOR USE: Polishing full or partial dentures

COLORS AND SIZES AVAILABLE:

ITEM #	GRIT	SIZE	PACKAGING
2011001S	FLOUR	1LB	JAR
2011001	FLOUR	25X1LB	JAR
2011005S	FLOUR	5LB	JAR
2011005	FLOUR	4X5LB	JAR
2011020	FLOUR	20LB	CARTON
2011040	FLOUR	40LB	CARTON
2013001S	FINE	1LB	JAR
2013001	FINE	25X1LB	JAR
2013005S	FINE	5LB	JAR
2013005	FINE	4X5LB	JAR
2013020	FINE	20LB	CARTON
2013040	FINE	40LB	CARTON
2015001S	MEDIUM	1LB	JAR
2015001	MEDIUM	25X1LB	JAR
2015005S	MEDIUM	5LB	JAR
2015005	MEDIUM	4X5LB	JAR
2015020	MEDIUM	20LB	CARTON
2015040	MEDIUM	40LB	CARTON
2016001S	COARSE	1LB	JAR
2016001	COARSE	25X1LB	JAR
2016005S	COARSE	5LB	JAR
2016005	COARSE	4X5LB	JAR
2016020	COARSE	20LB	CARTON
2016040	COARSE	40LB	CARTON
2017001S	X COARSE	1LB	JAR
2017001	X COARSE	25X1LB	JAR
2017005S	X COARSE	5LB	JAR
2017005	X COARSE	4X5LB	JAR
2017020	X COARSE	20LB	CARTON
2017040	X COARSE	40LB	CARTON
2018001S	XX COARSE	1LB	JAR
2018001	XX COARSE	25X1LB	JAR
2018005S	XX COARSE	5LB	JAR
2018005	XX COARSE	4X5LB	JAR
2018020	XX COARSE	20LB	CARTON
2018040	XX COARSE	40LB	CARTON



ABRASIVES

GLASS BEADS

Glass Beads are used to divest chrome cobalt, non-precious crown and bridge alloys, and all-ceramic cases. May also be used to polish chrome cobalt frameworks. Precision screened and desiccant packaged to prevent moisture.

INDICATION FOR USE: Divesting, ceramic cases, polishing chrome cobalt

COLORS AND SIZES AVAILABLE:

ITEM #	GRIT	SIZE	PACKAGING
2302025	50 MICRON	25LB	CARTON
2352025	50 MICRON	25LB	PAIL
2302050	50 MICRON	50LB	CARTON
2352050	50 MICRON	50LB	PAIL
2303525	80 MICRON	25LB	CARTON
2353525	80 MICRON	25LB	PAIL
2303550	80 MICRON	50LB	CARTON
2353550	80 MICRON	50LB	PAIL



WALNUT SHELLS

Walnut Shells from Garreco are ideal for use after divesting full and partial dentures to remove gypsum from the acrylic surface. 20/30 grit only.

INDICATION FOR USE: Removing gypsum from acrylic surfaces

COLORS AND SIZES AVAILABLE:

ITEM #	SIZE	PACKAGING
2400020	20LB	CARTON
2400040	40LB	CARTON





DYNASTY™ HI-IMPACT ACRYLIC

DYNASTY™ HI-IMPACT
Acrylic Denture Base



Dynasty Hi-Impact is a cadmium free denture base acrylic with ultimate strength and flexibility. It is dimensionally accurate, color stable, requires no change in technique, and has excellent working characteristics. Dynasty Hi-Impact polishes to a dense, high gloss finish.

INDICATION FOR USE: Denture base

CURE TYPE: Heat cure

INJECTABLE: Yes

MICROWAVABLE: Yes

LIQUID TYPE: Liquid Monomer #1

POWDER/LIQUID RATIO: 3:1

COLORS AND SIZES AVAILABLE:

PACKAGING	SHADE 99
1LB JAR	3748801S
4X1LB JAR	3748801
5LB JAR	3748805S
4X5LB JAR	3748805
20LB PAIL	3748820



DYNASTY™ HI-IMPACT ACRYLIC

DYNASTY™ HI-IMPACT POUR TYPE

Acrylic Denture Base

SELF
CURE

Dynasty Hi-Impact Pour Type is the first dental acrylic to achieve high strength and high impact resistance from a self cure formulation, achieving a fracture toughness 35% greater than conventional self cure denture base polymers. This new self cure material allows dental labs to cope with high production while keeping energy costs and labor to a minimum by avoiding the equipment, increased labor, energy use and cure times associated with heat cured materials.

INDICATION FOR USE: Denture base

CURE TYPE: Self cure

INJECTABLE: Yes

MICROWAVABLE: No

LIQUID TYPE: Liquid Monomer for Hi-Impact Pour Type Acrylics Only

POWDER/LIQUID RATIO: 2:1

COLORS AND SIZES AVAILABLE:

PACKAGING	SHADE 99
1LB JAR	3938801S
4X1LB JAR	3938801
5LB JAR	3938805S
4X5LB JAR	3938805
20LB PAIL	3938820



DYNASTY™ HI-IMPACT POUR TYPE LIQUID MONOMER

Dynasty Hi-Impact Pour Type Liquid Monomer has been specifically design to pair with the Dynasty Hi-Impact Pour Type Denture Base powder acrylic. This liquid formula is “non-yellowing” and packaged in ultraviolet proof containers to ensure the longest possible shelf life.

INDICATION FOR USE: Liquid for acrylic denture base powders

POWDER/LIQUID RATIO: 2:1

SIZES AVAILABLE:

SIZE	ITEM #
4OZ	3940004S
4X4OZ	3940004
8OZ	3940008S
4X8OZ	3940008
1QT	3940030S
4X1QT	3940030



AXCENT ACRYLIC™

PREMIUM Acrylic Denture Base

Premium Denture Base Acrylic may be compared with the best in the industry. This material offers long working time, strength, and flexibility for considerably less cost than competitors.

INDICATION FOR USE: Denture base

CURE TYPE: Heat cure

INJECTABLE: Yes

MICROWAVABLE: Yes

LIQUID TYPE: Liquid Monomer #1

POWDER/LIQUID RATIO: 3:1

COLORS AND SIZES AVAILABLE:



PACKAGING	PINK	LIGHT FIBER	DARK FIBER	LIGHT CHAR.	PURPLE FIBER	SHADE 99	CLEAR
1LB JAR	3108001S	3100201S	3100301S	3100401S	3103401S	3108801S	3100101S
4X1LB JAR	3108001	3100201	3100301	3100401	3103401	3108801	3100101
5LB JAR	3108005S	3100205S	3100305S	3100405S	3103405S	3108805S	3100105S
4X5LB JAR	3108005	3100205	3100305	3100405	3103405	2108805	3100105
20LB PAIL	3108020	3100220	3100320	3100420	3103420	3108820	3100120

STANDARD Acrylic Denture Base

A new manufacturing process enables Garreco to offer this impact resistant material at substantial cost savings over lower quality competitors.

INDICATION FOR USE: Denture base

CURE TYPE: Heat cure

INJECTABLE: Yes

MICROWAVABLE: Yes

LIQUID TYPE: Liquid Monomer #1

POWDER/LIQUID RATIO: 3:1

COLORS AND SIZES AVAILABLE:



PACKAGING	PINK	LIGHT FIBER	DARK FIBER	LIGHT CHAR.	PURPLE FIBER	SHADE 99	CLEAR
1LB JAR	3708001S	3700201S	3700301S	3700401S	3703401S	3708801S	3700101S
4X1LB JAR	3708001	3700201	3700301	3700401	3703401	3708801	3700101
5LB JAR	3708805S	3700205S	3700305S	3700405S	3703405S	3708805S	3700105S
4X5LB JAR	3708805	3700205	3700305	3700405	3703405	3708805	3700105
20LB PAIL	3708820	3700220	3700320	3700420	3703420	3708820	3700120

AXCENT ACRYLIC™

**15 MINUTE
Acrylic Denture Base**

**HEAT
CURE**

15-Minute Denture Base is a revolutionary acrylic material. Simply pack as usual, place into boiling water for 15 minutes, then finish as usual. Cures without porosity and its strength is comparable to most premium products.

INDICATION FOR USE: Denture base

CURE TYPE: 15 minute boil

INJECTABLE: No

MICROWAVABLE: No

LIQUID TYPE: Liquid Monomer #2

POWDER/LIQUID RATIO: 3:1

COLORS AND SIZES AVAILABLE:



PACKAGING	PINK	LIGHT FIBER	DARK FIBER	LIGHT CHAR.	PURPLE FIBER	SHADE 99	CLEAR
1LB JAR	3708001S	3730201S	3700301S	3730401S	3733401S	3738801S	3730101S
4X1LB JAR	3708001	3730201	3700301	3730401	3733401	3738801	3730101
5LB JAR	3708805S	3730205S	3700305S	3730405S	3733405S	3738805S	3730105S
4X5LB JAR	3708805	3730205	3700305	3730405	3733405	3738805	3730105
20LB PAIL	3708820	3730220	3700320	3730420	3733420	3738820	3730120

**REPAIR
Acrylic Denture Base**

**HEAT
CURE**

Repair Acrylic is color matched to all Garreco denture base materials. Contains special additives to prevent “yellowing” and may be used for any acrylic repair.

INDICATION FOR USE: Denture base repair

CURE TYPE: Heat cure

INJECTABLE: Yes

MICROWAVABLE: Yes

LIQUID TYPE: Liquid Monomer #1

POWDER/LIQUID RATIO: 3:1

COLORS AND SIZES AVAILABLE:



PACKAGING	PINK	LIGHT FIBER	DARK FIBER	LIGHT CHAR.	PURPLE FIBER	SHADE 99	CLEAR
1LB JAR	3018001S	3010201S	3010301S	3010401S	3013401S	3018801S	3010101S
4X1LB JAR	3018001	3010201	3010301	3010401	3013401	3018801	3010101
5LB JAR	3018005S	3010205S	3010305S	3010405S	3013405S	3018805S	3010105S
4X5LB JAR	3018005	3010205	3010305	3010405	3013405	3018805	3010105
20LB PAIL	3018020	3010220	3010320	3010420	3013420	3018820	3010120

AXCENT ACRYLIC™

POUR TYPE
Acrylic Denture Base



Recently reformulated, our Pour Type Acrylic has been specifically designed for use with all pour techniques. Ample pouring time and superior flexural strength are added benefits of this product.

INDICATION FOR USE: Denture base

CURE TYPE: Self cure

INJECTABLE: No

MICROWAVABLE: No

LIQUID TYPE: Liquid Monomer #5 (Pour Type)

POWDER/LIQUID RATIO: 2:1

COLORS AND SIZES AVAILABLE:



PACKAGING	LIGHT FIBER	SHADE 909
1LB JAR	3750201S	3750901S
4X1LB JAR	3750201	3750901
5LB JAR	3750205S	3750905S
4X5LB JAR	3750205	3750905
20LB PAIL	3750220	3750920

ORTHODONTIC
Acrylic Powder



Axcent Orthodontic Acrylic has been formulated with added stackability for orthodontic applications. Nine stock shades may be mixed to make hundreds of unique shades. Micro-fine metallic glitter may be added to any acrylic appliance for added characterization.

INDICATION FOR USE: Orthodontistry procedures

CURE TYPE: Self cure

INJECTABLE: No

MICROWAVABLE: No

LIQUID TYPE: Clear Powder - Liquid Monomer #4,

All Other Powders - Liquid Monomer #3

POWDER/LIQUID RATIO: 2:1

COLORS AND SIZES AVAILABLE:



PACKAGING	CLEAR	WHITE	BLACK
4OZ BOTTLE	3760104S	3763004S	3766304S
4X4OZ BOTTLE	3760104	3763004	3766304
1LB JAR	3760101S	3763001S	3766301S
4X1LB JAR	3760101	3763001	3766301
5LB JAR	3760105S	3763005S	3766305S
4X5LB JAR	3760105	3763005	3766305
20LB PAIL	3760120	3763020	3766320

AXCENT ACRYLIC™

TRAY MATERIAL

Acrylic Designed for Custom Trays



Tray Material by Garreco is auto curing and designed for custom trays.

INDICATION FOR USE: Custom trays

CURE TYPE: Self cure

INJECTABLE: No

MICROWAVABLE: No

LIQUID TYPE: Liquid Monomer #6 (Tray Liquid)

POWDER/LIQUID RATIO: 2:1

COLORS AND SIZES AVAILABLE:



PACKAGING	BLUE	WHITE
1LB JAR	3301501S	3303001S
4X1LB JAR	3301501	3303001
5LB JAR	3301505S	3303005S
4X5LB JAR	3301505	3303005
20LB PAIL	3301520	3303020

SOFTLINER RELINE MATERIAL KIT

Axcent™ Acrylic Softliner Reline Material is a high-quality permanent soft reline material. Placed in a new or old denture, it helps to improve the health of gum and tissue by absorbing some of the pressure of mastication (acts as a tissue-conditioning material) and helps to provide the maximum retention possible by utilizing undercuts in the bone and gum which hard liners may not be able to provide without causing irritation.

The Softliner Reline Base comes in a 2-part cartridge with mixing tips that negates manual measuring and mixing making easy and consistent results. The cartridge base also allows material to never be exposed to air, thus keeping unused product fresh and ready for next application. Long term expiration is appealing to all levels of production.

INDICATION FOR USE: Axcent™ Softliner is a silicone based liner to provide cushioning between the mouth and denture bases

KITS AVAILABLE:



ITEM #	ITEM DESCRIPTION
8601010	SOFTLINER RELINE MATERIAL KIT
8601030	SOFTLINER APPLICATION TOOLS KIT
8601020	SOFTLINER REFILL KIT
8601040	SOFTLINER PRIMER KIT
8601060	SOFTLINER GUN KIT
8601050	SOFTLINER VARNISH KIT

AXCENT ACRYLIC™

LIQUID MONOMER

Garreco Liquid Monomers have been specifically matched to optimally perform with all Garreco powders. All liquid formulas are cross linked, “non-yellowing”, and packaged in ultraviolet proof containers to ensure the longest possible shelf life.

INDICATION FOR USE: Liquid for acrylic powders

COLORS AND SIZES AVAILABLE:

ITEM #	TYPE	SIZE	PACKAGING
3780008S	HEAT CURE	8OZ	BOTTLE
3780008	HEAT CURE	4X8OZ	BOTTLE
3780030S	HEAT CURE	1QT	BOTTLE
3780030	HEAT CURE	4X1QT	BOTTLE
3800008S	15 MINUTE	8OZ	BOTTLE
3800008	15 MINUTE	4X8OZ	BOTTLE
3800030S	15 MINUTE	1QT	BOTTLE
3800030	15 MINUTE	4X1QT	BOTTLE
3830008S	SELF CURE	8OZ	BOTTLE
3830008	SELF CURE	4X8OZ	BOTTLE
3830030S	SELF CURE	1QT	BOTTLE
3830030	SELF CURE	4X1QT	BOTTLE
3860008S	HI-CLARITY	8OZ	BOTTLE
3860008	HI-CLARITY	4X8OZ	BOTTLE
3860030S	HI-CLARITY	1QT	BOTTLE
3860030	HI-CLARITY	4X1QT	BOTTLE
3910008	POUR TYPE	8OZ	BOTTLE
3910008S	POUR TYPE	4X8OZ	BOTTLE
3910030	POUR TYPE	1QT	BOTTLE
3910030S	POUR TYPE	4X1QT	BOTTLE
3920008	TRAY LIQUID	8OZ	BOTTLE
3920008S	TRAY LIQUID	4X8OZ	BOTTLE
3920030	TRAY LIQUID	1QT	BOTTLE
3920030S	TRAY LIQUID	4X1QT	BOTTLE
3845008S	SELF CURE - PINK	8OZ	BOTTLE
3845008	SELF CURE - PINK	4X8OZ	BOTTLE
3845030S	SELF CURE - PINK	1QT	BOTTLE
3845030	SELF CURE - PINK	4X1QT	BOTTLE



WHICH MONOMER DO I USE?

Use the chart below to determine which monomer to use with each acrylic powder.

ACRYLIC POWDER	LIQUID MONOMER
DYNASTY HI-IMPACT	HEAT CURE
DYNASTY HI-IMPACT POUR TYPE	HI-IMPACT POUR TYPE
PREMIUM	HEAT CURE
STANDARD	HEAT CURE
15 MINUTE	15 MINUTE
REPAIR	SELF CURE
POUR TYPE	POUR TYPE
ORTHODONTIC (CLEAR)	HI-CLARITY
ORTHODONTIC (COLORED)	SELF CURE
TRAY MATERIAL	TRAY LIQUID

FIND A GARRECO DEALER NEAR YOU

As a manufacturer of premium dental laboratory products, we are proud of our strong dealership network that spans multiple countries across the globe. Based on your location, one of our sales experts can provide you with a complete list of authorized Garreco dealers/distributors in your area.

To find an authorized Garreco Dealer near you, please call **1-800-334-1443 x806** or email us at sales@garreco.com





CROWN & BRIDGE

INVESTMENT


ULTRAVEST™ Investment Material

Ultravest excels when utilized with any crown & bridge alloy (except titanium) or any ceramic system. This phosphate bonded, graphite-free formula has a micro-fine particle matrix which allows precise attention to detail. This unique product may be used with either conventional or fast casting techniques and is designed for use with ring or ringless systems. UltraVest may be the only investment you will ever need!

INDICATION FOR USE: Investment for all crown & bridge alloys

ITEM #	SIZE	PACKAGING
6100012	50x100G	POUCHES
6109023	100x90G	POUCHES



 **WARNING:** Cancer - www.P65Warnings.ca.gov


PERFECT PRESS™ Investment Material

Perfect Press is a graphite-free, phosphate bonded precision investment developed especially for pressed ceramics. Utilizing a fast casting technique, Perfect Press not only divests easily, but delivers accurate consistent fits. Conventional or fast casting. For use with Hi-X or Regular investment/expansion liquids.

INDICATION FOR USE: Investment for all crown & bridge alloys

ITEM #	SIZE	PACKAGING
6150016	50x100G	POUCHES
6150026	100x100G	POUCHES



 **WARNING:** Cancer - www.P65Warnings.ca.gov

INVESTMENT

INVESTMENT LIQUID

For Ultravest™, Perfect Press™, and Vitex™ 3D Vest

Investment Liquid for use with Ultravest, Perfect Press, and Vitex™ 3D Vest. Available in 1 liter bottles in various quantities.

INDICATION FOR USE: Investment material liquid



ITEM #	SIZE	PACKAGING
6200003S	1L	BOTTLE
6200412	4x1L	BOTTLE
6201647	16x1L	BOTTLE

ULTRAVEST LIQUID REQUIREMENTS:	
Case Size	Liquid
50x100g	1Liter
100x90g	2Liters
150x60g	2Liters

PERFECT PRESS LIQUID REQUIREMENTS:	
Case Size	Liquid
50x100g	1Liter
100x100g	2Liters

3D VEST LIQUID REQUIREMENTS:	
Case Size	Liquid
50x400g	1Liter

FLEXI-RING™

Ringless Crown & Bridge System

Advanced Ringless System allows for complete expansion of investment materials due to flexible silicone construction. Each ring is constructed with inner ribs to provide additional surface area for optimum thermal conductivity. All rings are extremely durable and designed for years of continued use in all conditions. Available in 60g, 90g, 180g, and 300g sizes.

INDICATION FOR USE: Ringless investment procedures

ITEM #	DESCRIPTION	SIZE
6240001	EX-SMALL RING (2.5CM/1")	60G
6240101	SMALL RING (4.3CM/1 3/4")	90G
6240201	MEDIUM RING (6CM/2 3/8")	180G
6240301	LARGE RING (7.5CM/3")	300G
6260001	EX-SMALL BASE	60G
6260101	SMALL BASE	90G
6260201	MEDIUM BASE	180G
6260301	LARGE BASE	300G
6220005	STARTER KIT (INCLUDES 90G, 180G, AND 300G RINGS AND BASES)	ALL



C&B SUPPLIES

DIE SPACER

Precision Fit Agent

Die Spacer is used to provide sufficient space between the die and crown to accommodate cements used during seating. Garreco fast drying Die Spacers are formulated with additional ingredients to help seal and protect dies from wear and abrasion. Garreco's quality guarantees uniform thickness time after time. Available in 2 oz (brush cap) bottles or 4 oz (refill) bottles.

INDICATION FOR USE: Provides sufficient space between die and crown to accommodate cement

ITEM #	TYPE	SIZE	PACKAGING
7002002	BLUE SPACER	2OZ	BRUSH CAP BOTTLE
7002012	BLUE SPACER	12X2OZ	BRUSH CAP BOTTLE
7002004	BLUE SPACER	4OZ	REFILL BOTTLE
7002014	BLUE SPACER	12X4OZ	REFILL BOTTLE
7002102	GREY SPACER	2OZ	BRUSH CAP BOTTLE
7002112	GREY SPACER	12X2OZ	BRUSH CAP BOTTLE
7002104	GREY SPACER	4OZ	REFILL BOTTLE
7002114	GREY SPACER	12X4OZ	REFILL BOTTLE
7002202	THINNER	2OZ	BRUSH CAP BOTTLE
7002212	THINNER	12X2OZ	BRUSH CAP BOTTLE
7002204	THINNER	4OZ	REFILL BOTTLE
7002214	THINNER	12X4OZ	REFILL BOTTLE



HAZMAT ITEM:
Special shipping conditions may apply to these items.



WARNING: Cancer - www.P65Warnings.ca.gov

DEBUBBLIZER

Surfactant

A superior alcohol-free surfactant for reducing surface tension on wax patterns, impressions, and models.

INDICATION FOR USE: Reducing surface tension on wax patterns, impressions, and models

ITEM #	COLOR	SIZE	PACKAGING
3551508S	BLUE	8OZ	BOTTLE
3551508	BLUE	4x8OZ	BOTTLE
3551530S	BLUE	1QT	BOTTLE
3551530	BLUE	4x1QT	BOTTLE



PLASTER SEPARATOR

Apply as a paint-on or dip for film-free separation of plaster and stone.

INDICATION FOR USE: Plaster separator

ITEM #	COLOR	SIZE	PACKAGING
3488330S	RED	1QT	BOTTLE
3488330	RED	4x1QT	BOTTLE
3488332S	RED	1GAL	BOTTLE
3488332	RED	4x1GAL	BOTTLE



C&B SUPPLIES

VEST-IT™**Refractory Material**

Vest-it is an asbestos-free, high heat refractory material used to create custom firing pegs.

INDICATION FOR USE: Creating custom firing pegs

ITEM #	QUANTITIES
6840002	1 SYRINGE (10CC)

VOLUME DISCOUNT ON:

100 - 249 SYRINGES

250 - 499 SYRINGES

500+ SYRINGES

**HONEYCOMB TRAYS**

Specifically designed high heat firing trays for use with all furnaces and techniques. Available with ceramic or metal pins.

INDICATION FOR USE: Platform for firing crowns and bridges

ITEM #	STYLE	PACKAGING
7760801	CERAMIC PINS	2 TRAYS W/ PINS
7760802	METAL PINS	2 TRAYS W/ PINS

**FLEX TRAY**

This non-stick firing tray combines the best features of rigid trays and soft liners. You'll love the flexibility, durability, and multiple uses for this tray


INDICATION FOR USE: Flexible tray material for firing crowns and bridges

ITEM #	QUANTITIES
6840001	1 BAG (4 TRAYS/BAG)

VOLUME DISCOUNT ON:

25 - 99 BAGS

100+ BAGS

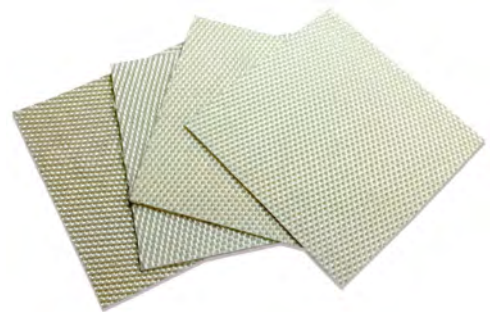**DISPOSABLE PLUNGERS**

Universal disposable plungers for pressing ceramics. Pre-heat not required. Available in 2 sizes (2 gram and 3 gram).

INDICATION FOR USE: Pressing ceramics

ITEM #	SIZE	PACKAGING
6080000	2 GRAM	50 PER CARTON
6080300	3 GRAM	50 PER CARTON



C&B SUPPLIES

CRUCIBLES

Type A, non-contaminating, quartz crucibles. Available in slotted or hooded. Sold in pairs.

INDICATION FOR USE: Melting alloys

ITEM #	TYPE & QUANTITIES
CRU-H	HOODED, 2/PACK
CRU-S	SLOTTED, 2/PACK



NATURAL AGATE

Agate Palette and Spatulas for Dental Porcelain

Polished agates provide an ideal platform for mixing dental porcelain. Their varicolor bands add beauty and uniqueness to these one of a kind natural palettes. Agates provide a calming influence, improving perception and concentration, thus enhancing the ceramist's creativity. Agates offer years of use because they do not abrade like porcelain or glass trays/slabs.

The Agate Spatula is an agate instrument ideal for mixing porcelain. The spatulas are metal-free to prevent metal from contaminating the porcelain, plus both ends can be used for mixing making it the perfect accompaniment to a natural agate palette.

INDICATION FOR USE: Mixing tray for dental porcelain

ITEM #	DESCRIPTION
6880010	AGATE PALETTE/SPATULA COMBO
680008	AGATE MIXING PALETTE
6850009	AGATE SPATULA



HUNDREDS OF AGATES TO CHOOSE FROM!

To view our Natural Agate Gallery, visit us at Garreco.com/products/crown-bridge/agate/ or scan this QR code now

PHYX-IT

Model and Die Adhesive

Phyx-It was developed to repair voids and undercuts in models and dies. Additional uses include: connecting sprues, connectors for bridges, adhering dowel pins, and attachment of plastic articulators.

INDICATION FOR USE: Repairs voids and undercuts in models and dies

ITEM #	DESCRIPTION
6910033	MODEL & DIE ADHESIVE (1/3 OZ)
6910002	SOLVENT (2 OZ)
6910102	ACCELERATOR (2 OZ)
6910016	ACCELERATOR (16 OZ)





GENERAL USE

SEPARATING FILM

Garreco Separating Film is a superior alternative to tinfoil and rinses completely from teeth and metal frameworks to provide unparalleled separation from the acrylic and mold. This product is available in Regular and Concentrate.

INDICATION FOR USE: Separating film

REGULAR

ITEM #	COLOR	SIZE	PACKAGING
3458330S	RED	1QT	BOTTLE
3458330	RED	4X1QT	BOTTLE
3458340S	RED	1GAL	BOTTLE
3458340	RED	4X1GAL	BOTTLE
3450130S	CLEAR	1QT	BOTTLE
3450130	CLEAR	4X1QT	BOTTLE
3450140S	CLEAR	1GAL	BOTTLE
3450140	CLEAR	4X1GAL	BOTTLE

CONCENTRATE

ITEM #	COLOR	SIZE	PACKAGING
3458330S	RED	1GAL	BOTTLE
3458330	RED	4X1GAL	BOTTLE



⚠ WARNING: Cancer - www.P65Warnings.ca.gov

LOOKING FOR PLASTER SEPARATOR ?

Plaster Separator is located on page 27; if you need more assistance please contact us at 800-334-1443 ext. 806


PERFORM™ PRODUCTS

PERFORM™ TOOTH SIEVE

PERform Tooth Sieves are a must when using pour techniques. These compartmented plastic sieves are heat resistant and durable for years of service.

INDICATION FOR USE: Used for easy scalding of teeth

ITEM #	COLOR
4100130	BLUE

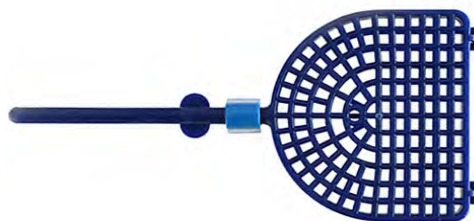

PERFORM™ TOOTH CUTTING BUR

PERform Tooth Cutting Burs are used for cutting diatorics in acrylic teeth.

INDICATION FOR USE: Cutting diatorics

ITEM #	SIZE
4100140	ONE SIZE



PERFORM™ FILTER PAPERS

PERform™ Filter Papers for use with PERform™ unit.

INDICATION FOR USE: Filters papers

ITEM #	QUANTITIES
4100100	500/PACK



CAN'T FIND WHAT YOU ARE LOOKING FOR?

For assistance with other PERform™ items or equipment, please contact us at 800-334-1443 ext. 806

SUBSCRIBE AND SAVE

Join our mailing list to get special sales and promotions delivered directly to your inbox. Visit Garreco.com/signup today to get started. Check out Garreco.com/promotions to see the current promotions we have available.



Garreco.com/signup



Garreco.com/promotions



PARTIAL DENTURES




VITEX™ PRODUCTS

VITEX™ INVESTMENT

Vitex (Ethyl Silicate) partial denture framework investment material has been developed especially for producing chrome cobalt castings. This product offers precise details due to an ultra fine particle size and unique two part blend of liquids. May be used for all high heat techniques. Formulated to mix and pour easily, Vitex investment maintains excellent integrity when used in conjunction with Vitex Model Hardener.

INDICATION FOR USE: Investment material

100LB KIT

ITEM #	DESCRIPTION	PACKAGING
5007000 	100LB KIT	3 CARTONS
Includes: 2 x 50lb Investment Sand, 2 x 118oz "A" Liquid, and 2 x 14.8oz "B" Liquid		
5012200 	100LB KIT LIQUID ONLY	CARTON
Includes: 2 x 118oz "A" Liquid, and 2 x 14.8oz "B" Liquid		
5008000	50LB SAND ONLY	CARTON
5009000 	118OZ "A" LIQUID	CARTON
5011000	14.8OZ "B" LIQUID	CARTON

50LB KIT


ITEM #	DESCRIPTION	PACKAGING
5006000	50LB KIT	2 CARTONS
Includes: 1 x 50lb Investment Sand, 4 x 28oz "A" Liquid, and 4 x 3.7oz "B" Liquid		
5010000	28OZ "A" LIQUID ONLY	CARTON
5012000	3.7OZ "B" LIQUID ONLY	CARTON



 Denotes an item that cannot be shipped by FedEx Ground; must be shipped FedEx Freight only



HAZMAT ITEM:
Special shipping conditions may apply to these items.

 **WARNING:** Cancer - www.P65Warnings.ca.gov


VITEX™ PRODUCTS

VITEX™ 3D VEST

Vitex™ 3D Vest is a phosphite bonded, graphite free and very fine grain precision investment especially for partial denture frameworks made with “Rapid-Prototyping” technologies. 3D Vest is specially developed to burn out plastic frames made with Rapid-Prototyping resins.

INDICATION FOR USE: Investment material designed for 3D Printed/Milled Wax and Plastic Scaffolding.

ITEM #	DESCRIPTION
5013400	50x400G POUCHES



 **WARNING:** Cancer - www.P65Warnings.ca.gov



VITEX™ REFRACTORY MODEL HARDENER

Vitex Refractory Model Hardener (“Cold Dip”) is designed to be used at room temperature with a warm model. Use of this material exhibits the same results as older techniques, while saving time and money.

INDICATION FOR USE: Model hardener

ITEM #	DESCRIPTION
5005030S	1QT BOTTLE
5005030	4X1QT BOTTLE
5005000 	1GAL BOTTLE
5005040 	4X1GAL BOTTLE

 Denotes an item that cannot be shipped by FedEx Ground; must be shipped FedEx Freight only



HAZMAT ITEM: Special shipping conditions may apply to these items.

VITEX™ DUPLICATING MATERIAL

Although developed for use in the construction of partial denture frameworks, this versatile hydrocolloid material has proven to be excellent for a variety of dental related tasks. With superior reproduction of detail and endurance for multiple melting cycles. Available in two types (all purpose and regular). All Purpose melts at 190° and holds at 135°. Regular melts at 200° and holds at 125°. Available in Neutral or Blue.

INDICATION FOR USE: Partial denture framework construction



ALL PURPOSE (Melt at 190°F, hold at 135°F)

ITEM #	COLOR	SIZE
5001015	BLUE	2 GALLON
5001050	BLUE	5 GALLON

REGULAR (Melt at 200°F, hold at 125°F)

ITEM #	COLOR	SIZE
5001515	BLUE	2 GALLON
5001550	BLUE	5 GALLON

VITEX™ PRODUCTS

VITEX™ ALLOY

This unique blend of metals provides short term flexibility and long term rigidity required by today's advanced technicians. Vitex partial denture alloy has been developed for use with all high heat techniques, and is priced well below other premium alloys.

INDICATION FOR USE: Dental alloy


































ITEM #	DESCRIPTION	PACKAGING
5018500	1LB BULLETS	CARTON



VITEX™ SUPERFORMS

Vitex Superform Patterns are made from a wax/plastic blend. This mix offers all the advantages of both components with optimum burnout. Available in 33 styles in protective packaging. 10 cards per box.

INDICATION FOR USE: Patterns for casting alloy

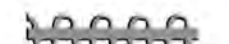 B-1 5051010 BOX	 B-2 5052010 BOX	 B-4 5054010 BOX	 B-6 5056010 BOX	 G-45 G-55 5083010 BOX 5084010 BOX
 GC-1 5085010 BOX	 RC-1 5103510 BOX	 GR-4 5099010 BOX	 RR-2 5107010 BOX	 FC-1 5081010 BOX
 GC-3 5087010 BOX	 GD-1 5093010 BOX	 GR-5 5101010 BOX	 BR-1 5072010 BOX	 BP-1 5068010 BOX
 GC-5 5089010 BOX	 GD-2 5094010 BOX	 GR-6 5102010 BOX	 BR-2 5073010 BOX	 BP-11 5071010 BOX
 GC-6 5091010 BOX	 GD-3 5095010 BOX	 GR-6 5102010 BOX	 RT-0 5109010 BOX	 BL-1 5063010 BOX
 GC-7 5092010 BOX	 GD-4 5097010 BOX	 CR-1 5076010 BOX	 RT-1 5111010 BOX	 BL-11 5065010 BOX
 GD-3L 5096010 BOX	 GD-5 5098010 BOX	 RI-1 5104010 BOX		



ROTARY

ROTARY

PERLADIA POLISHERS Diamond Impregnated and Long Lasting



Perlada Polishers have been the standard in porcelain polishing for many years and can be used with any porcelain. They are long lasting, generate no heat, and eliminate the need for re-glazing. No porcelain chipping and no polishing paste. All Perlada Polishers are impregnated with diamonds!

INDICATION FOR USE: Shaping, finishing and polishing with no chipping or heat generation.

25MM - KNIFE EDGE

	Ochre Extra Coarse Gross Trim Zirconia PDS25XCZO		Red Coarse Quick Reduction PDS25CR		Tan Medium- Coarse Contouring PDS25MCT		Blue Medium Finishing PDS25MB
	Green Fine Initial Polish PDS25FG		Purple Extra Fine Pre-Polish PDS25XFP		Yellow High Shine Glaze PDS25HSY		Black Coarse Zirconia PDS25CZBK

15MM - SQUARE EDGE

	Ochre Extra Coarse Gross Trim Zirconia PDS15XCZO		Red Coarse Quick Reduction PDS15CR		Grey Medium-Coarse Contouring PDS15MCGR		Blue Medium Finishing PDS15MB
---	--	---	---	---	--	---	--

ROTARY

NATURALGLAZE
Diamond Polishing System

Garreco ensures lasting good looks and the highest quality for all your restorations with Garreco rotary instruments. Naturalglaze polishers are excellent for shaping, finishing, and polishing ceramic materials with no chipping or heat generation.

INDICATION FOR USE: Shaping, finishing, and polishing ceramic with no chipping or heat generation.

25MM - KNIFE EDGE



17MM - SQUARE EDGE



FLAME POLISHERS



ROTARY - INTERPROXIMAL

Garreco has several dental rotary discs that are specifically designed to reach interproximal areas for Ceramics and Zirconia/Pressables.

INDICATION FOR USE: Shaping, finishing and polishing interproximal areas.

20MM - INTERPROXIMAL



VISIT GARRECO.COM TO VIEW AND DOWNLOAD THE COMPLETE ROTARY CATALOG



SILICONE PRODUCTS

**SILAFORM™
Silicone Putty**

Silaform silicone putty is a one to one mix material which is excellent for a variety of procedures. Silaform 90 (blue) has a Shore A hardness of 90 while Silaform 70 (purple) is more flexible and has a Shore A hardness of 70. Both materials resist deformation and abrasion normally associated with silicone putty. The absence of expansion and contraction offers endless possibilities for this economically priced product from Garreco.

INDICATION FOR USE: Duplicating material

MIXING RATIO: 1:1, A/B

SILAFORM™ 90 EXTRA HARD (BLUE)

ITEM #	SIZE	PACKAGING
6801505	1.5KG SET	CANISTER
6801555	25KG SET	PAIL

SILAFORM™ 70 HARD (PURPLE)

ITEM #	SIZE	PACKAGING
6803404	1.5KG SET	CANISTER
6803455	25KG SET	PAIL



SILICONE PRODUCTS

ACCUSIL™

Silicone Duplicating Material

A high quality platinum based silicone duplicating material for use with all techniques. Accuracy is constant and predictable at 0.0% expansion/contraction. Excellent strength of 16-18 Shore A hardness. May be poured multiple times with gypsum or investment without deterioration. Suitable for use in dispensers.

INDICATION FOR USE: Duplicating material

MIXING RATIO: 1:1, A/B

WORKING TIME: >5 min

CURING TIME: 30-45 min

ELONGATION OF TEAR: ~310%

ITEM #	SIZE	DESCRIPTION
6705005	1KG SET	1 ACTIVATOR & 1 BASE
6705055	25KG SET	1 ACTIVATOR & 1 BASE

STARTER KIT

ITEM #	DESCRIPTION
6700006	ACCUSIL STARTER KIT
Starter Kit Includes: 1 each 8oz Accusil A & B, 10" Duplicating Ring, Duplicating Cross, Silasavers (1 each: sm, med, lg)	



RAPIDSIL™

Silicone Duplicating Material

Lower cost alternative to Accusil with a slightly reduced Shore A hardness of 15 and a faster curing time. Rapidsil offers no expansion or contraction. Due to the use of a platinum base, this material may be poured hundreds of times with equally pleasing results.

INDICATION FOR USE: Duplicating material

MIXING RATIO: 1:1, A/B

WORKING TIME: >4 min

CURING TIME: 15-20 min

ELONGATION OF TEAR: ~300%

ITEM #	SIZE	DESCRIPTION
6701805	1KG SET	1 ACTIVATOR & 1 BASE
6701855	25KG SET	1 ACTIVATOR & 1 BASE



SILICONE PRODUCTS

SILSAVER™

Specially designed plastic spacers for silicone duplication. Adds stability to silicone, while saving material cost. Re-usable and available in 3 sizes.

INDICATION FOR USE: Spacers for silicone duplication



ITEM #	DESCRIPTION
6480101	SMALL SILSAVER (EACH)
6480201	MEDIUM SILSAVER (EACH)
6480301	LARGE SILSAVER (EACH)

DUPLICATING CROSS Placement Device

Re-usable, threaded placement device for SilsaVERS.

INDICATION FOR USE: Silsaver placement



ITEM #	DESCRIPTION
6440003	DUPLICATING CROSS

SPRUE FORMERS

Specially designed plastic sprue formers for use with refractory casting techniques.

INDICATION FOR USE: Refractory casting



ITEM #	DESCRIPTION
6600001	SPRUE FORMERS (10/PKG)

DUPLICATING RING

Disposable, clear, self adhesive duplicating ring, for use with silicone duplication techniques.

INDICATION FOR USE: Duplication techniques



ITEM #	DESCRIPTION
6460002	DUPLICATING RING (40M/ROLL)

SILICONE PRODUCTS

25 KG TOTE STAND

This custom-made adjustable tote stand is designed to hold 2 x 25kg silicone totes for use with our silicone dispenser. Hoses and fittings are included.

INDICATION FOR USE: Stand to hold 2 x 25kg silicone totes



ITEM #	DESCRIPTION
6720301	25KG TOTE STAND

SILICONE DISPENSER

Our silicone dispenser is the latest in European design and function. Simple one button mixing and dispensing makes this the easiest dispenser yet. Allows continuous dispensing of ready-mixed silicone by means of disposable mixing tips (unit comes with 100 mixing tips). Works with either 1kg bottles or 25kg totes.

INDICATION FOR USE: Continuous dispensing of silicone

ITEM #	DESCRIPTION
6720300	SILICONE DISPENSER (110V)



DISPOSABLE MIXING TIPS AVAILABLE TO ORDER FOR SILICONE DISPENSER

For more information regarding disposable mixing tips, please contact us at 800-334-1443 ext. 806



SAVE 59% OF YOUR SILICONE WITH OUR EASY MODEL DUPLICATING TECHNIQUE

We've created a step by step tutorial that shows you how to easily duplicate dental models and save 59% of your silicone using by Accusil/Rapidsil silicone, Garreco Duplicating Ring, Duplicating Cross, Silsaver, and gypsum. All of these products work together for an easy model duplication system that you can implement to save time and money in the lab.



Scan this QR code with your smartphone or visit <http://bit.ly/2zyj6EU> to view this technique.

HOW TO ORDER



FIND A GARRECO DEALER NEAR YOU

As a manufacturer of premium dental laboratory products, we are proud of our strong dealership network that spans multiple countries across the globe. Based on your location, one of our sales experts can provide you with a complete list of authorized Garreco dealers/distributors in your area.

To find a Garreco dealer near you, please call **1-800-334-1443 x806** to speak with a Garreco Sales Associate. You may also email us at **orders@garreco.com** and we will promptly reply within 1 business day (usually in 1 hr or less during business hours).

View our website at **www.garreco.com** for our complete line of dental stones and plasters, acrylics, abrasives, rotary, and dental lab supplies.

CONTACT US

Normal Operating Hours:

Monday-Friday, from 8:00am to 5:00pm

Physical Address:

Garreco, LLC
430 Hiram Road
Heber Springs, AR 72543

Mailing Address:

Garreco, LLC
PO Box 1258
Heber Springs, AR 72543

Toll Free Number: (800) 334-1443

Local Number: (501) 362-6261

Fax Number: (501) 362-2264

Online:

Website: Garreco.com

Blog: Garreco.com/blog

Email Addresses:

General Info: info@garreco.com

Sales: sales@garreco.com

Support: support@garreco.com

Orders: orders@garreco.com


Social Media:

Facebook: Facebook.com/GarrecoLLC/

Twitter: Twitter.com/GarrecoDental

LinkedIn: LinkedIn.com/company/garreco-dental

REFERRAL PROGRAM

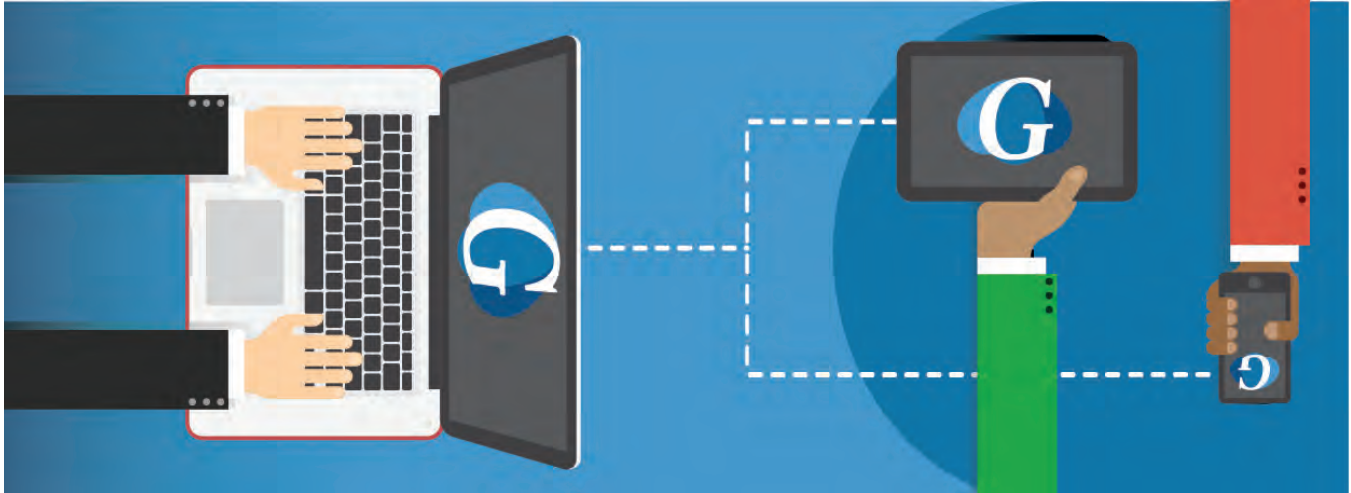

REFER A LAB PROGRAM

Refer a Dental Lab to Garreco and you each get \$150!*


Your endorsement matters to us. When a lab requests your advice in finding new products, Garreco Dental appreciates your recommendation. Refer another Dental Lab to Garreco and we'll reward you both with \$150 each if your referral invests \$1,500 or more in Garreco products (within 3 months).

TO PARTICIPATE, FOLLOW THESE SIMPLE STEPS:

1. Find a Dental Lab to refer to Garreco.
2. Provide and fill out the information on the form located at Garreco.com/refer.
3. Send the digital referral form to us or via email, fax or postal mail.
4. Start ordering Garreco products through your dealers.
5. Receive \$150 for you and your referral when they invest a total of \$1,500 in Garreco products in a 3 month period.

VISIT GARRECO.COM/REFER TODAY TO GET STARTED NOW!

Valid from 8/01/2017. Referral's \$1,500 net sales investment must be achieved within 3 months of submitting the referral. Send digital referral form to Garreco Dental, PO Box 1258, Heber Springs, Arkansas 72543 or fax to 501-362-6261 or email sales@garreco.com. These deals cannot be combined with any other promotion or discount. Other restrictions may apply. Garreco reserves the right to change or cancel this program at any time. Not responsible for typographical errors. See website for more details.

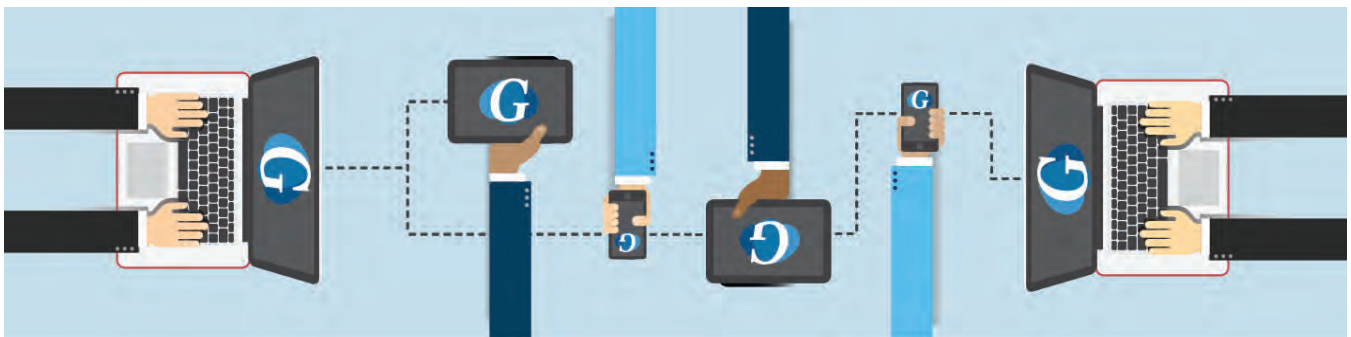

VISIT US ONLINE AT WWW.GARRECO.COM

SAFETY DATA SHEETS

Safety Data Sheets
**All files in Adobe Acrobat PDF format*

ABRASIVES

- Aluminum Oxide, Brown
- Aluminum Oxide, White
- Glass Beads
- Pumice
- Resilience
- Walnut Shells

ACRYLICS

- Accent Acrylics 15 Minute Denture Base Powder
- Accent Acrylics Dynasty Denture Base Powder
- Accent Acrylics Liquid Monomer #1, Heat Cure
- Accent Acrylics Liquid Monomer #2, HC 15-Minute
- Accent Acrylics Liquid Monomer #3, Self Cure
- Accent Acrylics Liquid Monomer #4, HC Hi-Clearity

SUPPORT
Go to Garreco.com/support to view/download our Safety Data Sheets

PROMOTIONS

SILAFORM SILICONE PUTTY

BUY 2, GET 1 FREE ON SILAFORM SET (A & B)

Silafom silicone putty is a one to one mix material which is excellent for a variety of procedures. Silafom 90 (blue) has a Shore A hardness of 90 while Silafom 70 (purple) is more flexible and has a Shore A hardness of 70. Both materials resist deformation and abrasion normally associated with silicone putty. The absence of expansion and contraction offers endless possibilities for this economically priced product from Garreco.

[SEE OFFER](#)

CURRENT PROMOTIONS

2nd Quarter Promotions
EFFECTIVE FROM APRIL 1 - JUNE 30, 2017
*Some restrictions apply

GARRECO PROMOTIONS
Go to Garreco.com/promotions to view our current promotions and special deals.

GARRECO DENTAL BLOG

To get the latest news from Garreco, check out Garreco.com/blog and join in on the conversation!

LET'S KEEP IN TOUCH

At Garreco.com/signup you can subscribe to our mailing list to get all of the latest from Garreco Dental delivered to your inbox.



Use this QR Code with your smartphone to go to the signup page now!



PRODUCT INDEX

#	
15 Minute Denture Base	21
3D Vest Investment Material	33
A	
ABRASIVES	14
Accusil™, Silicone Duplicating Material	38
ACRYLICS	18
Agate, Natural	29
Alloy, Vitex™	34
Aluminum Oxide	15
Apex™, Type IV-R Resin Reinforced Die Stone	4
Articulating Plaster, Type II Gypsum	12
Axcent™ Acrylic	20
B	
Base Stones	8
C	
C&B Supplies	27
Cad-Scan™, Type IV-R Resin Die Stone	5
CROWN & BRIDGE	25
Crucibles, Quartz	29
Custom Products	3
D	
Debubblizer	27
Denture Products	32
Diamond Impregnated Polishers	35
Die Spacer	27
Die Stones	4
Dispenser, Silicone	40
Disposable Plungers	28
Duplicating Cross	39
Duplicating Material, Vitex™	33
Duplicating Products	37
Duplicating Ring	39
Dynasty™ Hi-Impact Denture Base	18
Dynasty™ Hi-Impact Pour Type Denture Base	19
E	
Excalibur™, Type IV High Strength Die Stone	5
Expansion Stone	11
Express Set™, Type IV Specialty Stone	6
F	
Fastone™, Type III Base Stone	9
Filter Papers, PERform™	31
Flask It™, Type III Base Stone	10
Flasking Stone™, Type III Base Stone	10
Flex Tray™	28
Flexi-Ring™ Ringless Crown & Bridge System	26
G	
General Use	30
Glass Beads	17
GYP SUM	4
H	
Hi-Impact Acrylic	18
Hi-Impact Pour Type Liquid Monomer	19
Honeycomb Trays	28
I	
Interproximal Rotary	36
Investment Liquid	26
Investment Material, Vitex™ 3D Vest	33
Investment Material, Vitex™	32
Investment Products	25
L	
LAB SUPPLIES	30
Labstone™, Type III Base Stone	9
Liquid Monomer	24
M	
Model and Die Adhesive	29
Mounting Stone™, Type III Base Stone	10
N	
NaturalGlaze™ Diamond Polishing System	36
NaturalRock™ EH, Type IV-R Resin Reinforced	6
O	
Operator Gypsum	13
Orthodontic Acrylic Powder	22
Orthodontic Gypsum	11
Orthopaster™, Type III Base Stone	11
Orthostone™, Type III Base Stone	11
P	
PARTIAL DENTURES	32
PD Polishers	35
Perfect Press™, Investment Material	25
PERform™ Products	31
Perladia Polishers	35
Phyx-It Model and Die Adhesive	29
Plaster	12
Plaster Separator	27
Plaster, Type II Gypsum	12
Pour Type Denture Base	22
Premium Denture Base	20
Pristine® II, Type II Gypsum	13
Pristine® III, Type III Base Stone	13
Pristine® V, Type IV Die Stone	13
Product Samples	3
Pumice	16
R	
Rapidsil™, Silicone Duplicating Material	38
Refractory Material, Vest-It™	28
Refractory Model Hardener, Vitex™	33
Reline Material Kit	23
Repair Denture Base	21
Resilience™, Polishing and Finishing Media	14
ROTARY	35
S	
Sample, Product	3
Satinstone™, Type IV Die Stone	6
Separating Film	30
Silaform™	37
SILICONE	37
Silicone Dispenser	40
Silicone Duplicating Material	38
Silicone Putty	37
Silsaver™	39
Softliner Reline Material Kit	23
Specialty Gypsum	11
Sprue Formers	39
Standard Denture Base	20
Superforms, Vitex™	34
T	
Tecstone™ FL - Flowable, Type III Base Stone	8
Tecstone™, Type III Base Stone	8
Tooth Cutting Bur, PERform™	31
Tooth Sieve, PERform™	31
Tote Stand, 25KG	40
Tray Material	23
U	
Ultravest™, Investment Material	25
V	
Veri-Die™ 20, Type IV High Expansion Die Stone	7
Veri-Die™, Type V Die Stone	7
Vest-It™, Refractory Material	28
Vitex™ Products	32
W	
Walnut Shells	17
X	
X-Roc™	11



CONTACT US



Normal Operating Hours:
Monday-Friday, 8:00am to 5:00pm



Physical Address:
Garreco, LLC
430 Hiram Road
Heber Springs, AR 72543



Mailing Address:
Garreco, LLC
PO Box 1258
Heber Springs, AR 72543



Toll Free Number: (800) 334-1443
Local Number: (501) 362-6261
Fax Number: (501) 362-2264



Online:
Website: Garreco.com
Blog: Garreco.com/blog



Email Addresses:
General Info: info@garreco.com
Sales: sales@garreco.com
Support: support@garreco.com
Orders: orders@garreco.com



Social Media:
Facebook: Facebook.com/GarrecoLLC/



Twitter: Twitter.com/GarrecoDental