



# Garreco

*Quality Dental Lab Supplies Since 1988*

3D RESINS   GYPSUM   ABRASIVES   ACRYLICS   LAB SUPPLIES   ROTARY   WAX



## PRODUCT CATALOG

PAGE  
**04**



**3D RESINS**

PAGE  
**06**



**DENTAL STONE & GYPSUM**

PAGE  
**12**



**ABRASIVES**

PAGE  
**14**



**ACRYLICS**

PAGE  
**17**



**CROWN & BRIDGE**

PAGE  
**19**



**LAB SUPPLIES**

PAGE  
**20**



**PARTIAL DENTURES**

PAGE  
**23**



**ROTARY**

PAGE  
**26**



**SILICONES**

PAGE  
**28**



**WAXES**

## INTRODUCTION & RESOURCES

About Garreco .....	1
Find a Dealer .....	1
Custom Products .....	3
Product Samples .....	3
Online Resources .....	3

## 3D RESINS

Model 3D Pro Beige .....	4
Model 3D Pro Grey .....	4
Model 3D Ultra .....	4
Gingiva 3D .....	5
WaxCast 3D .....	5
Clean99Isopropyl Alcohol .....	5

## GYPSUM

<b>Die Stones</b> .....	6
Apex™, Type IV-R Resin .....	6
Cad-Scan™, Type IV-R Resin .....	6
Excalibur™, Type IV High Strength .....	7
Express Set™, Type IV Specialty .....	7
NaturalRock EH™, Type IV-R Resin .....	7
Satinstone™, Type IV .....	7
Veri-Die™, Type IV .....	8
Veri-Die™ 20, Type V High Expansion .....	8
<b>Base Stones</b> .....	8
Tecstone™, Type III .....	8
Tecstone™ FL - Flowable, Type III .....	8
Fastone™, Type III .....	9
Labstone™, Type III .....	9
Flask It™, Type III .....	9
Flasking Stone™, Type III .....	9
Mounting Stone™, Type III .....	9
<b>Orthodontic Gypsum</b> .....	10
Orthopaster™, Type III .....	10
Orthostone™, Type III .....	10
<b>Specialty Gypsum</b> .....	10
X-Roc™, Expansion Stone .....	10
<b>Plaster</b> .....	10
Plaster, Type II .....	11
Articulating Plaster, Type II .....	11
<b>Operatory Gypsum</b> .....	11
Pristine® II, Type II .....	11
Pristine® III, Type III .....	11
Pristine® V, Type IV .....	11

## ABRASIVES

Aluminum Oxide .....	12
Glass Beads .....	12
Pumice .....	13
Resilience™ Acrylic Polishing Media .....	13
Walnut Shells .....	13

## ACRYLICS

<b>Dynasty™ Hi-Impact Acrylic</b> .....	14
Hi-Impact Denture Base .....	14
<b>Axcent™ Acrylic</b> .....	14
Premium Denture Base .....	14
Standard Denture Base .....	14
15 Minute Denture Base .....	15
Repair Denture Base .....	15
Pour Type Denture Base .....	15
Orthodontic Acrylic Powder .....	16
Tray Material .....	16
Liquid Monomer .....	16

## CROWN & BRIDGE

<b>Investment Products</b> .....	17
Ultravest™, Investment Material .....	17
Perfect Press™, Investment Material .....	17
Investment Liquid .....	17
<b>C&amp;B Supplies</b> .....	17
Crucibles .....	17
Die Spacer .....	18
Disposable Plungers .....	18
Flex Tray .....	18
Debubbler .....	19
Vest-It™ Refractory Material .....	19
Honeycomb Trays .....	19
Flexi-Ring™ Ringless C&B System .....	19

## LAB SUPPLIES

Ultrasonic Cleaner #1 .....	19
Ultrasonic Cleaner #2 .....	19
Ultrasonic Cleaner #5 .....	19
Separating Film .....	20
<b>PERform™ Products</b> .....	20
Tooth Sieve .....	20
Tooth Cutting Bur .....	20
Filter Papers .....	20

## PARTIAL DENTURES

<b>Vitex™ Products</b> .....	20
Investment Material .....	20
3D Vest Investment Material .....	21
Refractory Model Hardener .....	21
Duplicating Material .....	21
Alloy .....	21
Superforms .....	22

## ROTARY

Perlada (PD) Polishers .....	23
NaturalGlaze Diamond Polishing System .....	23
Carbides - FG & HP .....	24
Hub Brushes .....	24
Muslin Buffs .....	25
Bristle Brushes .....	25
Interproximal Rotary .....	26

## SILICONE

Silaform™ Silicone Putty .....	26
Accusil™ Silicone Duplicating Material .....	26
Rapidsil™ Silicone Duplicating Material .....	27
Silsaver™ .....	27
Duplicating Cross .....	27
Sprue Formers .....	27
Duplicating Ring .....	27

## DENTAL WAX

<b>Base Plate Waxes</b> .....	28
Medium Soft Base Plate Wax .....	28
Medium Hard Base Plate Wax .....	28
All Season Base Plate Wax .....	28
<b>Orthodontic Wax</b> .....	28
Ortho Wax Compacts .....	28
<b>Utility Wax</b> .....	29
Regular Boxing Wax .....	29
Wax Square Ropes .....	29
Large Wax Utility Strips .....	29
Small Wax Utility Strips .....	29
Periphery Wax Strips .....	29

# INTRODUCTION

## ABOUT GARRECO

Garreco Dental has been providing quality dental lab products since 1988. As one of the leading manufacturers of dental gypsum and abrasive products, we take great pride in the quality and consistency of all of our products.

Our manufacturing process is ISO 13485 certified and uses only the best raw materials which we carefully formulate for exacting product standards and superior performance. In addition, 100% of our gypsum products are tested prior to being shipped to ensure quality and consistency.

In addition to our core dental business, Garreco also produces products for the fields of jewelry casting, microdermal abrasion, forensics, taxidermy, special effects, and hobby. Please contact our sales team for more information regarding these products.

Our dedicated team of professionals is committed to providing superior customer service and on-time delivery for all our products. We encourage you to reach out and allow us to impress you with our products, honest approach to business, and the Garreco Advantage.



## THE GARRECO ADVANTAGE

- ★ Quality Dental Lab Products since 1988
- ★ ISO Certified Manufacturer
- ★ Highest Quality Product for a Great Price
- ★ Consistency in Every Box
- ★ Right Product Delivered on Time
- ★ Custom Blending & Packaging

## HOW TO ORDER



## FIND A GARRECO DEALER NEAR YOU

As a manufacturer of premium dental laboratory products, we are proud of our strong dealership network that spans multiple countries across the globe. Based on your location, one of our sales experts can provide you with a complete list of authorized Garreco dealers/distributors in your area.

To find a Garreco dealer near you, please call **1-800-334-1443** to speak with a Garreco Sales Associate. You may also email us at **orders@garreco.com** and we will promptly reply within 1 business day (usually in 1 hr or less during business hours).



## SALES AND RETURN POLICIES

### TERMS

Our normal terms are Net-30 days. However, sales will be handled on a COD/Prepay basis until an accepted credit application is on file. See Pricing for additional information.

### MINIMUM ORDER

Our minimum order amount is \$100.00 for combined products per invoice, excluding freight charges. A \$5.00 handling fee will be assessed on all orders below stated minimum.

### PRICING

We make every attempt to reduce costs and keep our prices as low as possible. However, we find that incremental increases are necessary. Several factors are taken into account such as: raw material and packaging cost, transportation cost, volume per item, volume per order, and general cost of doing business.

### FREIGHT CHARGES

Because freight costs are often a large percentage of your total order, we endeavor to obtain the lowest shipping cost possible. Orders will be shipped, freight prepaid and all charges will be added to invoice, unless prior arrangements have been made.

### PALLET FEE

A pallet fee will be charged on all palletized orders.

### DAMAGE CLAIMS

Any damaged goods should be immediately noted on delivery receipt when freight is received (for FedEx deliveries, please notify driver). If damage is not noted, Garreco will not be responsible for corrective action. Damage claims must be submitted to Garreco within 7 days of delivery. Furthermore, Garreco will not be responsible for corrective action for carriers other than those preferred by Garreco.

### PRODUCT RETURNS

Returns on custom products, colors, and packaging are not allowed. Returns for stocked products are acceptable up to one year from invoice date with a RA number. Credit is subject to a 15% restock fee for Garreco Brand and 20% restock fee for Private Label. All returns received in unsatisfactory condition (damaged) will be rejected.



## CUSTOM PRODUCTS

See what Garreco Dental can do for you! As a private label manufacturer for some of the largest dental companies in the world, Garreco is well-versed in offering a premium product customized to your specific needs. We offer an in-house design team to design your private label artwork.



### CUSTOM COLOR

Our Custom Color Program offers a choice from a wide variety of pastel shades for their gypsum products, including matching company colors.



### CUSTOM SET

Custom Set Products are available for almost all Garreco gypsums. This service offers a specified set time from 2 to 40 minutes. Whether you need to finish one case quickly or a hundred cases at the same time, we can customize a product for you.



### CUSTOM PACKAGING & LABEL DESIGN

Custom Packaging and Private Label Design services are available for all Garreco products.


## PRODUCT SAMPLES

Want to try a Garreco Dental product before ordering? We offer samples for most of our products available upon request.

To make your sample request or get more information about our samples or products, please contact the Garreco Sales Team at **800-334-1443** or **email us at [sales@garreco.com](mailto:sales@garreco.com)**.



## ONLINE RESOURCES

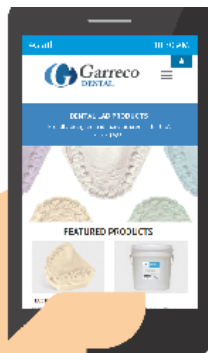

### GARRECO.COM & GARRECOPRINT.COM

Visit our websites to view our Featured Products and get technical support, download sales flyers and safety data sheets, and more. We have Safety Data Sheets available to download for each product. Our new website, [garrecoprint.com](http://garrecoprint.com) has everything you need to get started with our line of 3D resins.

### GARRECO DENTAL BLOG

Join us at [Garreco.com/blog](http://Garreco.com/blog) to read about what's happening at Garreco and get involved in the discussion!

### STAY IN TOUCH ON SOCIAL

Visit us on the social channels below:

**Facebook** - [Facebook.com/GarrecoMfg](https://Facebook.com/GarrecoMfg)

**Instagram** - [Instagram.com/garrecodental](https://Instagram.com/garrecodental)

**YouTube** - [Youtube.com/channel/UCHZtuMaiGnf6OEzgC8TXICQ](https://Youtube.com/channel/UCHZtuMaiGnf6OEzgC8TXICQ)

**Twitter/X** - [Twitter.com/GarrecoDental](https://Twitter.com/GarrecoDental)

**LinkedIn** - [Linkedin.com/company/garreco-dental](https://Linkedin.com/company/garreco-dental)

# GARRECO PRINT



## Save Money and Time with Garreco Print

Our new digital line is perfect for 3D printing a range of dental models including ortho, aligner, C&B and implant models, plus soft tissue and intricate parts. Garreco Print resins are open source and compatible on a wide range of popular 3D printers and UV-curing units.

Visit [Garrecoprint.com](http://Garrecoprint.com) to or scan the QR code to learn more about how you can save money and time with Garreco Print.



## 3D RESINS

### MODEL 3D PRO BEIGE & MODEL 3D PRO GREY - Photopolymer 3D Resin

Model 3D Pro Beige is a matte, beige opaque colored photopolymer. Its matte surface finish accentuates depth and detail in dental models, which accommodates a perfect view on undercuts and other important markers. In addition, 3D printed models from Model 3D Pro are excellent for crown & bridge build-ups. 3D-printed parts from this material have exceptional dimensional stability and extremely low shrinkage during printing. Model 3D Pro is easy to use on all open source SLA, MSLA and DLP 3D printers in the range of 385 – 420nm. This material has excellent properties like low shrinkage and low odor, accuracy and dimensional stability, making it perfect for the production of dental models, C&B, thermoforming and aligner models.

#### 3D-PRINTER COMPATIBILITY

Asiga Max and 4K, Ackuretta Dentiq, Photon and Phrozen Series, Open 385 - 420nm LCD /MSLA, DLP and SLA 3D printers. Visit [garrecoprint.com](http://garrecoprint.com) for the full list of compatible printers and settings, plus compatible cure boxes and post-processing instructions.

#### COMPATIBLE CURE BOXES

Ackuretta Curie, Formlabs Form Wash & Cure, Anycubic Wash & Cure 2.0, Asiga Flash Cure, Phrozen Cure Luna, and more. See our complete list of compatible Cure Boxes.

#### MODEL 3D PRO:

Beige 1kg (1000g) Bottle - 6851002S

Grey 1kg (1000g) Bottle - 6854502S



### MODEL 3D ULTRA - Photopolymer 3D Resin

Model 3D Ultra is an opaque peach colored photopolymer resin. Its color accentuates depth in your dental models, which accommodates a perfect view on undercuts in 3D printed models. 3D printed parts from this material have great dimensional stability and low shrinkage during printing. Model 3D Ultra is easy to use on all open source LCD and DLP 3D printers in the range of 385 – 420nm. This material has excellent properties like low shrinkage, low odor and high print speed making it perfect for the production of orthodontic and aligner models.

#### 3D-PRINTER COMPATIBILITY

Asiga Max and 4K, Shining3d AccuFab-D1, Photon and Phrozen Series, All low power, open source 385 -420nm LCD and DLP 3D-printers. Visit [garrecoprint.com](http://garrecoprint.com) for the full list of compatible printers and settings, plus compatible cure boxes and post-processing instructions.

#### COMPATIBLE CURE BOXES

Ackuretta Curie, Formlabs Form Wash & Cure, Anycubic Wash & Cure 2.0, Asiga Flash Cure, Phrozen Cure Luna, and more. See our complete list of compatible Cure Boxes.

#### MODEL 3D ULTRA:

1kg (1000g) Bottle - 6865502S



## 3D RESINS

### GINGIVA 3D - Photopolymer 3D Resin

Gingiva 3D is a pink colored photopolymer resin. It is a soft and exible material that can be used in combination with Model 3D Pro Beige or Grey material. Gingiva 3D opens the possibility to create the soft part for highly precise models, bridges and implant superstructures. Gingiva 3D is easy to use on open MSLA / LCD, DLP and SLA 3D printers in the range of 385 - 420nm. With its soft and exible properties, it is perfect for the production of gum parts and soft tissue which supports dental implant models and C&B projects.

#### 3D-PRINTER COMPATIBILITY

Asiga Series, Photon Series, Phrozen Series, Open 385 - 420nm LCD / MSLA, DLP and SLA 3D printers. Visit [garrecoprint.com](http://garrecoprint.com) for the full list of compatible printers and settings, plus compatible cure boxes and post-processing instructions.

#### COMPATIBLE CURE BOXES

Akuretta Curie, Formlabs Form Wash & Cure, Anycubic Wash & Cure 2.0, Asiga Flash Cure, Phrozen Cure Luna, and more. See our complete list of compatible Cure Boxes.

#### GINGIVA 3D:

1kg (1000g) Bottle - 6885002S



### WAXCAST 3D - Photopolymer 3D Resin

WaxCast 3D is a blue wax-based photopolymer with reliable processability and accuracy on LCD/MSLA, DLP and laser based 3D printers. 3D printed parts from this material captures intricate features and show crisp details. This wax-based material oers smooth surfaces with clean burnout for a reliable casting process. Create custom made elegant organic geometries according to your preferences with positive and negative engravings. WaxCast 3D creates perfect casting patterns for jewelry, dental and industrial parts.

#### 3D-PRINTER COMPATIBILITY

Asiga Series, Phrozen Sonic Mini 4K, Anycubic Photon, Elegoo Mars, Open 385 - 420nm LCD and SLA 3D-printers. Visit [garrecoprint.com](http://garrecoprint.com) for the full list of compatible printers and settings, plus compatible cure boxes and post-processing instructions.

#### COMPATIBLE CURE BOXES

Akuretta Curie, Formlabs Form Wash & Cure, Anycubic Wash & Cure 2.0, Asiga Flash Cure, Phrozen Cure Luna, and more. See our complete list of compatible Cure Boxes.

#### WAXCAST 3D:

1kg (1000g) Bottle - 6871502S



## POST-PROCESSING & CLEANING

### CLEAN 99 - Pure Isopropyl Alcohol

Clean 99 is a high-quality, 99% pure isopropyl alcohol for post-processing and cleaning 3D printed materials as well as general use purposes. This pure alcohol will eectively remove support structures, residue, and contaminants, leaving your prints clean and ready. Say goodbye to the hassle of cleaning and post-processing your 3D printed materials. Our Clean 99 Alcohol is the perfect companion for your printing needs, ensuring a seamless and easy experience every time.

**INDICATION FOR USE:** For cleaning and post-processing 3D printed material and general purposes uses.

#### CLEAN 99 PURE ISOPROPYL ALCOHOL:

1qt Bottle - 6790030S





# DENTAL STONES & PLASTERS

## Quality American Made Dental Stones & Plasters

We manufacture a full line of high-quality dental gypsums and plasters in a variety of popular set times, expansions, strength, and colors.

Products with custom colors, set times, expansions, and packaging are available. Please contact sales@garreco.com for more information on custom products.



## DIE STONES

### APEX™ RZ 10 & RZ 15 - Type IV-R Resin Reinforced Die Stone

Apex is a resin reinforced natural gypsum die material, available in two expansion variations (RZ 10, RZ 15) making it the perfect choice to complement any technique. Perfect for all crown and bridge applications as well as implant or precision attachment cases. Resin additives eliminate brittleness and add superior external hardness. The exclusive Apex flow agent allows for excellent stackability, while promoting free material flow during vibration.



**INDICATION FOR USE:** Crown and bridge; porcelain; implants; cast partials; high-precision cases.

**SETTING TIME:** 10-13 min

**WATER/POWDER RATIO:** 19ml:100g

**COMPRESSIVE STRENGTH:** 18,500 psi

**EXPANSION:** 0.10% (RZ 10), 0.15% (RZ 15)

#### RZ 10 COLORS AND SIZES:

Blue 25lb Box - 1301525 Blue 47.5lb Box - 1301550	Ivory 25lb Box - 1303525 Ivory 47.5lb Box - 1303550
Green 25lb Box - 1302525 Green 47.5lb Box - 1302550	White 25lb Box - 1303025 White 47.5lb Box - 1303050
Pink 25lb Box - 1305025 Pink 47.5lb Box - 1305050	Tan 25lb Box - 1306525 Tan 47.5lb Box - 1306550

#### RZ 15 COLORS AND SIZES:

Blue 25lb Box - 1321525 Blue 47.5lb Box - 1321550	Ivory 25lb Box - 1323525 Ivory 47.5lb Box - 1323550
Green 25lb Box - 1322525 Green 47.5lb Box - 1322550	White 25lb Box - 1323025 White 47.5lb Box - 1323050
Grey 25lb Box - 1324525 Grey 47.5lb Box - 1324550	Tan 25lb Box - 1326525 Tan 47.5lb Box - 1326550

 **WARNING:** Cancer - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)

### CAD-SCAN™ - Type IV-R Resin Die Stone

CAD-Scan is a unique gypsum material that contains a special scanning additive and resin reinforcement. CAD-Scan reflects light nicely for essential readability, allowing for a shorter working time, and is manufactured from an ultra-fine gypsum base, which prevents nicks and abrasions caused by brittleness. CAD-Scan has a low expansion and high compressive strength, making it the material of choice for all CAD/CAM restorations. It is perfect for production environments, as well as those single cases requiring a rapid set.

**INDICATION FOR USE:** Crown and bridge; porcelain; implants; cast partials; high-precision cases

**SETTING TIME:** 11-13 min

**WATER/POWDER RATIO:** 20ml:100g

**COMPRESSIVE STRENGTH:** 18,500 psi

**EXPANSION:** 0.06%-0.09%

#### CAD-SCAN COLORS AND SIZES:

Golden 10lb Canister - 1466710 Golden 25lb Box - 1466725 Golden 47.5lb Box - 1466750	Green 10lb Canister - 1462510 Green 25lb Box - 1462525 Green 47.5lb Box - 1462550	Buff 10lb Canister - 1461010 Buff 25lb Box - 1461025 Buff 47.5lb Box - 1461050
--	---	--



 **WARNING:** Cancer - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)

## DIE STONES

### EXCALIBUR™ - Type IV High Strength Die Stone

With its smooth mixing characteristics and ample setting time, Excalibur produces a “textbook” model every time. Die trimming becomes a worry-free procedure due to the ultra-fine gypsum base, which prevents nicks and abrasions caused by brittleness. Recommended for use as crown and bridge or partial denture models. Suitable for use with all impression materials. Gypsum hardeners are not recommended.

**INDICATION FOR USE:** Crown and bridge; porcelain; implants; cast partials; high-precision cases

**SETTING TIME:** 12-15 min

**WATER/POWDER RATIO:** 22ml:100g

**COMPRESSIVE STRENGTH:** 18,000 psi

**EXPANSION:** 0.09%

#### COLORS AND SIZES:

Blue 25lb Box - 1351525  
Blue 47.5lb Box - 1351550

Violet 25lb Box - 1352025  
Violet 47.5lb Box - 1352050

Tan 25lb Box - 1356525  
Tan 47.5lb Box - 1356550

Green 25lb Box - 1352525  
Green 47.5lb Box - 1352550

Grey 25lb Box - 1354525  
Grey 47.5lb Box - 1354550

White 25lb Box - 1353025  
White 47.5lb Box - 1353050



### EXPRESS SET™ - Type IV Specialty Stone

ExpressSet is a Type IV specialty stone with an extremely fast set time, smooth mixing consistency, and uniform color matrix. Models can be released from the impression in 5 minutes, producing a precision cast expressly suited for mouth guards, custom trays, study casts, models for opto-electronic/laser scanning with powder, temporaries and denture repairs.



**INDICATION FOR USE:** Crown and bridge; porcelain; implants; cast partials; high-precision cases

**SETTING TIME:** 60-90 seconds

**WATER/POWDER RATIO:** 23ml:100g

**COMPRESSIVE STRENGTH:** 14,000 psi

**EXPANSION:** 0.15%

#### COLORS AND SIZES:

Buff 25lb Box - 1401525  
Buff 47.5lb Box - 1401550

Pink 25lb Box - 1402525  
Pink 47.5lb Box - 1402550

### NATURALROCK™ EH - Type IV-R Resin Reinforced Die Stone

NaturalRock EH is a traditional, creamy, smooth, resin reinforced die stone which exhibits high compressive strength and very low expansion. This economical natural stone was designed to compare to Fujirock Die Stone and Resin Rock.



**INDICATION FOR USE:** Crown and bridge; porcelain; implants; cast partials; high-precision cases

**SETTING TIME:** 10-13 min

**WATER/POWDER RATIO:** 20ml:100g

**COMPRESSIVE STRENGTH:** 18,000 psi

**EXPANSION:** 0.08%

#### COLORS AND SIZES:

Light Tan 25lb Box - 1446725  
Light Tan 47.5lb Box - 1446750

### SATINSTONE™ - Type IV Die Stone

Comparable to Whip Mix's Silky-Rock, SatinStone™ is distinguished by exceptional toughness and an extraordinary “satin-smooth” finish that works well with all types of impression materials. Introducing only the finest raw materials allows this product to have a creamy and effortless mix each and every time. SatinStone™ also has a low expansion which makes it one of the most the most unique die stones on the market.



**INDICATION FOR USE:** Crown and bridge; porcelain; implants; cast partials; high-precision cases

**SETTING TIME:** ~10 min

**WATER/POWDER RATIO:** 23ml:100g

**COMPRESSIVE STRENGTH:** 15,000 psi

**EXPANSION:** 0.09%-0.11%

#### COLORS AND SIZES:

Rose 25lb Box - 1175025  
Rose 47.5lb Box - 1175050

White 25lb Box - 1173025  
White 47.5lb Box - 1173050

## DIE STONES

### VERI-DIE™ - Type IV Die Stone

Produced to handle the many different laboratory techniques that are being used in today's dental labs. With its range in application and versatility, Veri-Die is perfect for almost every requirement. Its physical properties allow Veri-Die to conform to the requirements for each use. Due to its range in particle size, Veri-Die virtually eliminates normal chipping and abrasion. Veri-Die has become the worry free material of choice by providing the technician ample working time and a moderate expansion rate.



**INDICATION FOR USE:** Crown and bridge; porcelain; implants; cast partials; high-precision cases

**SETTING TIME:** 10-13 min  
**WATER/POWDER RATIO:** 21ml:100g  
**COMPRESSIVE STRENGTH:** 16,500 psi  
**EXPANSION:** 0.15%

#### COLORS AND SIZES:

Blue 25lb Box - 1401525  
 Blue 47.5lb Box - 1401550

Green 25lb Box - 1402525  
 Green 47.5lb Box - 1402550

Pink 25lb Box - 1405025  
 Pink 47.5lb Box - 1405050

Ivory\* 25lb Box - 1403525  
 Ivory\* 47.5lb Box - 1403550

**WARNING:** Cancer - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)

### VERI-DIE™ 20 - Type V High Expansion Die Stone

Produced to handle the many different laboratory techniques that are being used in today's dental labs. With its range in application and versatility, Veri-Die is perfect for almost every requirement. Its physical properties allow Veri-Die to conform to the requirements for each use. Due to its range in particle size, Veri-Die virtually eliminates normal chipping and abrasion. Veri-Die has become the worry free material of choice by providing the technician ample working time and a moderate expansion rate.



**INDICATION FOR USE:** Crown and bridge; porcelain; implants; cast partials; high-precision cases

**SETTING TIME:** 10-13 min  
**WATER/POWDER RATIO:** 21ml:100g  
**COMPRESSIVE STRENGTH:** 18,000 psi  
**EXPANSION:** 0.20%

#### COLORS AND SIZES:

Blue 25lb Box - 1431225  
 Blue 47.5lb Box - 1431250

Green 25lb Box - 1432525  
 Green 47.5lb Box - 1432550

Cream 25lb Box - 1431425  
 Cream 47.5lb Box - 1431450

## BASE STONES

### TECSTONE™ - Type III Base Stone

Widely recognized as one of the best stones on the market and designed for ultimate versatility, Tecstone is recommended for use in models and flasking, in the production of full or partial dentures, and in basing crown and bridge models. A creamy mix, moderate compressive strength, and superior outward hardness are all characteristics of this unique formula.



**INDICATION FOR USE:** Full or partial denture models; flasking; basing crown and bridge models

**SETTING TIME:** 13-15 min  
**WATER/POWDER RATIO:** 30ml:100g  
**COMPRESSIVE STRENGTH:** 9,000 psi  
**EXPANSION:** 0.12%

#### COLORS AND SIZES:

Blue 25lb Box - 1011525  
 Blue 47.5lb Box - 1011550

Buff 25lb Box - 1011025  
 Buff 47.5lb Box - 1011050

White 25lb Box - 1013025  
 White 47.5lb Box - 1013050

### TECSTONE™ FL FLOWABLE - Type III Base Stone

This unique dental stone from Garreco is specifically designed for basing crown & bridge models. Tecstone FL – Flowable pours into base formers without vibration and yields a bubble-free base. Tecstone FL – Flowable far exceeds the strength of conventional base-stones while offering a dense, smooth result. Accuracy is also addressed by a low expansion.



**INDICATION FOR USE:** Basing crown and bridge models

**SETTING TIME:** 9-11 min  
**WATER/POWDER RATIO:** 23ml:100g  
**COMPRESSIVE STRENGTH:** 12,500 psi  
**EXPANSION:** 0.08%

#### COLORS AND SIZES:

Blue 25lb Box - 1031525  
 Blue 47.5lb Box - 1031550

Buff 25lb Box - 1031025  
 Buff 47.5lb Box - 1031050

White 25lb Box - 1033025  
 White 47.5lb Box - 1033050

**WARNING:** Cancer - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)

## BASE STONES

### FASTONE™ - Type III Base Stone




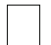
This dental stone offers the same unique characteristics as our Tecstone product, with a faster setting time. Used for full or partial denture models, denture flasking procedures, and recommended for use in basing crown and bridge models.



**INDICATION FOR USE:** Full or partial denture models; flasking; basing crown and bridge models

**SETTING TIME:** 10-12 min  
**WATER/POWDER RATIO:** 30ml:100g  
**COMPRESSIVE STRENGTH:** 9,000 psi  
**EXPANSION:** 0.12%

#### COLORS AND SIZES:

	Blue 25lb Box - 1101525 Blue 47.5lb Box - 1101550		Golden 25lb Box - 1106725 Golden 47.5lb Box - 1106750
	Buff 25lb Box - 1101025 Buff 47.5lb Box - 1101050		White 25lb Box - 1103025 White 47.5lb Box - 1103050

### LABSTONE™ - Type III Base Stone




The “work horse” of the gypsum family. Designed for ultimate versatility. Labstone is recommended for use in models and flasking, in the production of full or partial dentures, and in basing crown and bridge models. Creamy mix, moderate compressive strength, and superior outward hardness, are all characteristics of this unique formula.



**INDICATION FOR USE:** Full or partial denture models, denture flasking procedures, basing crown and bridge models

**SETTING TIME:** 11-14 min  
**WATER/POWDER RATIO:** 30ml:100g  
**COMPRESSIVE STRENGTH:** 8,500 psi  
**EXPANSION:** 0.12%

#### COLORS AND SIZES:

	Blue 25lb Box - 1051525 Blue 47.5lb Box - 1051550
	Buff 25lb Box - 1051025 Buff 47.5lb Box - 1051050
	White 25lb Box - 1053025 White 47.5lb Box - 1053050

### FLASK IT™ - Type III Base Stone


This base stone allows the technician to make the 2nd and 3rd pours during flasking without vibration. The result is an exceptionally smooth surface which greatly reduces air pockets, bubbles, and valuable time spent removing gypsum around teeth and under-cut areas. This is a must-try for all denture labs.



**INDICATION FOR USE:** Flasking procedures

**SETTING TIME:** 9-12 min  
**WATER/POWDER RATIO:** 33ml:100g  
**COMPRESSIVE STRENGTH:** 9,000 psi  
**EXPANSION:** 0.09%

#### COLORS AND SIZES:

	Pink 25lb Box - 1235025 Pink 47.5lb Box - 1235050
---	--

### FLASKING STONE™ - Type III Base Stone


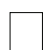
A moderate strength, low expansion stone, formulated for denture flasking procedures. Allows rapid heat transfer during boil-out, while maintaining optimum strength during packing. Divests quickly and easily.



**INDICATION FOR USE:** Flasking procedures

**SETTING TIME:** 9-12 min  
**WATER/POWDER RATIO:** 33ml:100g  
**COMPRESSIVE STRENGTH:** 9,000 psi  
**EXPANSION:** 0.09%

#### COLORS AND SIZES:

	Buff 25lb Box - 1201025 Buff 47.5lb Box - 1201050
	White 25lb Box - 1203025 White 47.5lb Box - 1203050

### MOUNTING STONE™ - Type III Base Stone



A high strength, fast setting, low expansion, articulating stone. Originally developed for implant cases where a low expansion articulating material is desired. Many professionals choose this product as a rapid setting alternative material for denture models and crown and bridge bases.



**INDICATION FOR USE:** Full or partial denture models, denture flasking procedures, basing crown and bridge models

**SETTING TIME:** 3-5 min  
**WATER/POWDER RATIO:** 30ml:100g  
**COMPRESSIVE STRENGTH:** 9,000 psi  
**EXPANSION:** 0.09%

#### COLORS AND SIZES:

	Grey 25lb Box - 1154525 Grey 47.5lb Box - 1154550
	White 25lb Box - 1153025 White 47.5lb Box - 1153050



## ORTHODONTIC STONE

### ORTHOPLASTER - Type III Base Stone

High strength plaster designed for orthodontic applications. Orthopaster combines smooth mixing, low expansion and ample setting time. Contains special whitening additives for a brilliant white color.



**INDICATION FOR USE:** Orthodontic applications

**SETTING TIME:** 10-13 min (Regular);  
6-9 (Fast)

**WATER/POWDER RATIO:** 35ml:100g

**COMPRESSIVE STRENGTH:** 6,500 psi

**EXPANSION:** 0.12%

**COLORS AND SIZES:**

☐ White Regular 25lb Box - 1600025  
White Regular 47.5lb Box - 1600050

☐ White Fast 25lb Box - 1620025  
White Fast 47.5lb Box - 1620050

 **WARNING:** Cancer - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)

### ORTHOSTONE - Type III Base Stone

OrthoStone exhibits the strength and accuracy of a quality base stone, plus special additives resulting in a brilliant white color. The superior color, working properties, and hardness make this the product of choice for those specializing in orthodontics.



**INDICATION FOR USE:** Orthodontic applications

**SETTING TIME:** 13-15 min

**WATER/POWDER RATIO:** 30ml:100g

**COMPRESSIVE STRENGTH:** 9,000 psi

**EXPANSION:** 0.12%

**COLORS AND SIZES:**

☐ White 25lb Box - 1253025  
White 47.5lb Box - 1253050

 **WARNING:** Cancer - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)

## SPECIALTY STONE

### X-ROC™ - Expansion Stone

X-Roc was developed for pouring the duplicate model in processing thermoplastic appliances such as: Flexite®, Valplast®, Flex-Star V, TCS®, FRS®, DuraFlex®, & DurAcetal®. This stone compensates for the shrinkage of these materials, thereby eliminating tight fits on removable prosthetics processed with conventional stones.



**INDICATION FOR USE:** Specialty duplication material

**SETTING TIME:** 22-25 min

**WATER/POWDER RATIO:** 29ml:100g

**COMPRESSIVE STRENGTH:** 13,000 psi

**EXPANSION:** 1.10%

**COLORS AND SIZES:**

☐ White 25lb Box - 1483025

## PLASTER

### PLASTER - Type II Gypsum

Dental Plasters offered by Garreco are formulated especially for use in today's complex field of dentistry. Designed for articulating and flasking procedures, Garreco plasters mix smoothly and exhibit excellent working characteristics. Available in three formulations with varying set times.



**INDICATION FOR USE:** Articulating and flasking procedures

**SETTING TIME:** 17-19 min (Slow);  
13-15 min (Regular); 9-11 min (Fast)

**WATER/POWDER RATIO:** 45ml:100g

**COMPRESSIVE STRENGTH:** 3,500 psi

**EXPANSION:** 0.19%

**PLASTER SLOW:**

25lb Carton - 1550025

47.5lb Carton - 1550050

**PLASTER FAST:**

25lb Carton - 1500025

47.5lb Carton - 1500050

**PLASTER REGULAR:**

25lb Carton - 1450025

47.5lb Carton - 1450050

## PLASTER

### ARTICULATING PLASTER - Type II Gypsum



Articulating Plaster exhibits similar characteristics to our Dental Plasters, but with a greatly reduced set time. Recommended for articulation use.

**INDICATION FOR USE:** Articulation

**SETTING TIME:** 3-5 min  
**WATER/POWDER RATIO:** 45ml:100g  
**COMPRESSIVE STRENGTH:** 3,500 psi  
**EXPANSION:** 0.19%

**COLORS AND SIZES:**

☐ White 25lb Box - 1560025  
☐ White 47.5lb Box - 1560050

## OPERATORY STONE

### PRISTINE® II - Type II Gypsum



A fast setting Type II model and articulating plaster exhibiting excellent working characteristics. Has a fresh, clean scent to mask odors picked up by impression materials which can be transferred to working models. Comes in stackable pails with flip top re-closable lid for easy use.

**INDICATION FOR USE:** Gypsum designed for dental office applications

**SETTING TIME:** 6-10 min  
**WATER/POWDER RATIO:** 45ml:100g  
**COMPRESSIVE STRENGTH:** 3,500 psi  
**EXPANSION:** 0.19%

**COLORS AND SIZES:**

☐ White 20lb Pail - 1500020

### PRISTINE® III - Type III Base Stone

With a broad range of applications this Type III stone may be used for full and partial denture models, denture flasking procedures, and for basing crown and bridge models. Has a fresh, clean scent to mask odors picked up by impression materials which can be transferred to working models. Comes in stackable pails with flip top re-closable lid for easy use.



**INDICATION FOR USE:** Gypsum designed for dental office applications

**SETTING TIME:** 8-10 min  
**WATER/POWDER RATIO:** 30ml:100g  
**COMPRESSIVE STRENGTH:** 9,000 psi  
**EXPANSION:** 0.12%

**COLORS AND SIZES:**

☒ Golden 25lb Pail - 1026725  
☐ White 25lb Pail - 1023025

### PRISTINE® V - Type IV Die Stone

Versatility for all your fixed and removable models. This Type IV die stone resists normal chipping and abrasion. Has a fresh, clean scent to mask odors picked up by impression materials which can be transferred to working models. Comes in stackable pails with flip top re-closable lid for easy use.



**INDICATION FOR USE:** Gypsum designed for dental office applications

**SETTING TIME:** 10-13 min  
**WATER/POWDER RATIO:** 21ml:100g  
**COMPRESSIVE STRENGTH:** 15,000 psi  
**EXPANSION:** 0.15%

**COLORS AND SIZES:**

☒ Green 25lb Pail - 1382525  
☐ Pink 25lb Pail - 1385025

# ABRASIVES

## ABRASIVES

### ALUMINUM OXIDE

Garreco Aluminum Oxide blast abrasives are among the finest in the world and used for divesting and cleaning of appliances. Our 99.6% pure (white) and 96% pure (brown) electro-fused Aluminum Oxide abrasives are machine-sized and computer-scanned to provide uniform particle size. Available in popular sizes and packed with desiccant to prevent moisture absorption.

**INDICATION FOR USE:** Divesting/cleaning of appliances

#### 25 MICRON WHITE:

25lb Carton - 2101025  
25lb Pail - 2151025  
47.5lb Carton - 2101050   
44.5lb Pail - 2151050

#### 60 MICRON WHITE:

25lb Carton - 2103025  
25lb Pail - 2103025  
47.5lb Carton - 2103050   
44.5lb Pail - 2153050

#### 250 MICRON WHITE:

25lb Carton - 2106025  
25lb Pail - 2156025  
47.5lb Carton - 2106050   
44.5lb Pail - 2156050

#### 50 MICRON WHITE:

2lb Pail - 2172002S  
6x2lb Pail - 2172002  
25lb Carton - 2102025  
25lb Pail - 2152025  
47.5lb Carton - 2102050   
44.5lb Pail - 2152050

#### 100 MICRON WHITE:

25lb Carton - 2104025  
25lb Pail - 2154025  
47.5lb Carton - 2104050   
44.5lb Pail - 2154050

#### 110 MICRON WHITE:

25lb Carton - 2104525  
25lb Pail - 2154525  
47.5lb Carton - 2104550   
44.5lb Pail - 2154550

#### 50 MICRON BROWN:

25lb Carton - 2202025  
25lb Pail - 2252025  
47.5lb Carton - 2202050   
44.5lb Pail - 2252050

 Denotes an item that cannot be shipped via FedEx Ground; must be shipped FedEx Freight only



### GLASS BEADS

Glass Beads are used to divest chrome cobalt, non-precious crown and bridge alloys, and all-ceramic cases. May also be used to polish chrome cobalt frameworks. Precision screened and desiccant packaged to prevent moisture.

**INDICATION FOR USE:** Divesting, ceramic cases, polishing chrome cobalt

#### 50 MICRON WHITE:

25lb Carton - 2302025  
25lb Pail - 2352025  
47.5lb Carton - 2302050  
44.5lb Pail - 2352050

#### 80 MICRON WHITE:

25lb Carton - 2303525  
25lb Pail - 2353525  
47.5lb Carton - 2303550  
44.5lb Pail - 2353550



## FIND A GARRECO DEALER NEAR YOU

As a manufacturer of premium dental laboratory products, we are proud of our strong dealership network that spans multiple countries across the globe. Based on your location, one of our sales experts can provide you with a complete list of authorized Garreco dealers in your area.

To find an authorized Garreco Dealer near you, please call **1-800-334-1443** or email us at [sales@garreco.com](mailto:sales@garreco.com)



## ABRASIVES

### PUMICE

A consistently fast, smooth cutting, long lasting abrasive. Suitable for polishing full or partial dentures. Available in a variety of particle sizes.

**INDICATION FOR USE:** Polishing full or partial dentures

#### FLOUR:

1lb Jar - 2011001S  
25x1lb Jar - 2011001  
5lb Jar - 2011005S  
4x5lb Jar - 2011005  
20lb Carton - 2011020  
40lb Carton - 2011040

#### FINE:

1lb Jar - 2013001S  
25x1lb Jar - 2013001  
5lb Jar - 2013005S  
4x5lb Jar - 2013005  
20lb Carton - 2013020  
40lb Carton - 2013040

#### MEDIUM:

1lb Jar - 2015001S  
25x1lb Jar - 2015001  
5lb Jar - 2015005S  
4x5lb Jar - 2015005  
20lb Carton - 2015020  
40lb Carton - 2015040

#### COARSE:

1lb Jar - 2016001S  
25x1lb Jar - 2016001  
5lb Jar - 2016005S  
4x5lb Jar - 2016005  
20lb Carton - 2016020  
40lb Carton - 2016040

#### X COARSE:

1lb Jar - 2017001S  
25x1lb Jar - 2017001  
5lb Jar - 2017005S  
4x5lb Jar - 2017005  
20lb Carton - 2017020  
40lb Carton - 2017040

#### XX COARSE:

1lb Jar - 2018001S  
25x1lb Jar - 2018001  
5lb Jar - 2018005S  
4x5lb Jar - 2018005  
20lb Carton - 2018020  
40lb Carton - 2018040



### RESILIENCE™ - Polishing/Finishing Media

Resilience™ is an all-purpose acrylic polishing media designed for faster and smoother polishing of acrylic appliances. Specially formulated with unique abrasives to cut faster, easier, and smoother than pumice-type products with no scratching or grooving. Resilience saves time and does not require a change in your current technique. Kind to your hands and does not give off a bad odor commonly associated with pumice.

**INDICATION FOR USE:** Polishing acrylic appliances

#### SIZES:

2.5lb Bottle - 2019003  
10lb Carton - 2019010  
20lb Carton - 2019020

 **WARNING:** Cancer - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)



#### BENEFITS OF USING RESILIENCE:

- It's easy to use
- Replacement for pumice media
- One grit size does it all
- No odor or bacteria
- No scratching or grooving
- Saves time
- High shine polishing
- Compares to Polycril and Dazzle

### WALNUT SHELLS

Walnut Shells from Garreco are ideal for use after divesting full and partial dentures to remove gypsum from the acrylic surface. 20/30 grit only.

**INDICATION FOR USE:** Removing gypsum from acrylic surfaces

#### SIZES:

20lb Carton - 2400020  
40lb Carton - 2400040





# ACRYLICS

## DYNASTY™ HI-IMPACT ACRYLIC

### DYNASTY™ HI-IMPACT - Acrylic Denture Base

Dynasty Hi-Impact is a cadmium free denture base acrylic with ultimate strength and flexibility. It is dimensionally accurate, color stable, requires no change in technique, and has excellent working characteristics. Dynasty Hi-Impact polishes to a dense, high gloss finish.

**INDICATION FOR USE:** Denture base

**CURE TYPE:** Heat cure

**INJECTABLE:** Yes

**MICROWAVABLE:** Yes

**LIQUID TYPE:** Monomer #1 Heat Cure

**POWDER/LIQUID RATIO:** 3:1

**SHADE 99:**

1lb Jar - 3748801S  
4x1lb Jar - 3748801  
5lb Jar - 3748805S  
4x5lb Jar - 3748805  
20lb Pail - 3748820



## AXCENT ACRYLIC™

### PREMIUM - Acrylic Denture Base

Premium Denture Base Acrylic may be compared with the best in the industry. This material offers long working time, strength, and flexibility for considerably less cost than competitors.

**INDICATION FOR USE:** Denture base

**CURE TYPE:**

Heat cure

**INJECTABLE:** Yes

**MICROWAVABLE:** Yes

**LIQUID TYPE:**

Monomer #1 Heat Cure

**POWDER/LIQUID**

**RATIO:** 3:1



**PINK:**

1lb Jar - 3108001S  
4x1lb Jar - 3108001  
5lb Jar - 3108005S  
4x5lb Jar - 3108005  
20lb Pail - 3108020

**LIGHT FIBER:**

1lb Jar - 3100201S  
4x1lb Jar - 3100201  
5lb Jar - 3100205S  
4x5lb Jar - 3100205  
20lb Pail - 3100220

**DARK FIBER:**

1lb Jar - 3100301S  
4x1lb Jar - 3100301  
5lb Jar - 3100305S  
4x5lb Jar - 3100305  
20lb Pail - 3100320

**PURPLE FIBER:**

1lb Jar - 3103401S  
4x1lb Jar - 3103401  
5lb Jar - 3103405S  
4x5lb Jar - 3103405  
20lb Pail - 3103420

**LIGHT CHARACTERIZED:**

1lb Jar - 3100401S  
4x1lb Jar - 3100401  
5lb Jar - 3100405S  
4x5lb Jar - 3100405  
20lb Pail - 3100420

**SHADE 99:**

1lb Jar - 3108801S  
4x1lb Jar - 3108801  
5lb Jar - 3108805S  
4x5lb Jar - 3108805  
20lb Pail - 3108820

**CLEAR:**

1lb Jar - 3100101S  
4x1lb Jar - 3100101  
5lb Jar - 3100105S  
4x5lb Jar - 3100105  
20lb Pail - 3100120

HEAT  
CURE

### STANDARD - Acrylic Denture Base

A new manufacturing process enables Garreco to offer this impact resistant material at substantial cost savings over lower quality competitors.

**INDICATION FOR USE:** Denture base

**CURE TYPE:**

Heat cure

**INJECTABLE:** Yes

**MICROWAVABLE:** Yes

**LIQUID TYPE:**

Monomer #1 Heat Cure

**POWDER/LIQUID**

**RATIO:** 3:1



**PINK:**

1lb Jar - 3708001S  
4x1lb Jar - 3708001  
5lb Jar - 3708805S  
4x5lb Jar - 3708805  
20lb Pail - 3708820

**LIGHT FIBER:**

1lb Jar - 3700201S  
4x1lb Jar - 3700201  
5lb Jar - 3700205S  
4x5lb Jar - 3708805  
20lb Pail - 3700220

**DARK FIBER:**

1lb Jar - 3700301S  
4x1lb Jar - 3700301  
5lb Jar - 3700305S  
4x5lb Jar - 3700305  
20lb Pail - 3700320

**PURPLE FIBER:**

1lb Jar - 3703401S  
4x1lb Jar - 3703401  
5lb Jar - 3703405S  
4x5lb Jar - 3703405  
20lb Pail - 3703420

**LIGHT CHARACTERIZED:**

1lb Jar - 3700401S  
4x1lb Jar - 3700401  
5lb Jar - 3700405S  
4x5lb Jar - 3700405  
20lb Pail - 3700420

**SHADE 99:**

1lb Jar - 3708801S  
4x1lb Jar - 3708801  
5lb Jar - 3708805S  
4x5lb Jar - 3708805  
20lb Pail - 3708820

**CLEAR:**

1lb Jar - 3700101S  
4x1lb Jar - 3700101  
5lb Jar - 3700105S  
4x5lb Jar - 3700105  
20lb Pail - 3700120

HEAT  
CURE

## AXCENT ACRYLIC™

### 15 MINUTE - Acrylic Denture Base

15-Minute Denture Base is a revolutionary acrylic material. Simply pack as usual, place into boiling water for 15 minutes, then finish as usual. Cures without porosity and its strength is comparable to most premium products.

**INDICATION FOR USE:**  
Denture base

**CURE TYPE:**  
15 minute boil  
**INJECTABLE:** No  
**MICROWAVABLE:** No  
**LIQUID TYPE:**  
Monomer #2 15 Minute  
**POWDER/LIQUID RATIO:** 3:1



#### PINK:

1lb Jar - 3738001S  
4x1lb Jar - 3738001  
5lb Jar - 3738805S  
4x5lb Jar - 3738805  
20lb Pail - 3738820

#### PURPLE FIBER:

1lb Jar - 3733401S  
4x1lb Jar - 3733401  
5lb Jar - 3733405S  
4x5lb Jar - 3733405  
20lb Pail - 3733420

#### LIGHT FIBER:

1lb Jar - 3730201S  
4x1lb Jar - 3730201  
5lb Jar - 3730205S  
4x5lb Jar - 3708805  
20lb Pail - 3730220

#### LIGHT CHARACTERIZED:

1lb Jar - 3730401S  
4x1lb Jar - 3730401  
5lb Jar - 3730405S  
4x5lb Jar - 3730405  
20lb Pail - 3730420

#### CLEAR:

1lb Jar - 3730101S  
4x1lb Jar - 3730101  
5lb Jar - 3730105S  
4x5lb Jar - 3730105  
20lb Pail - 3730120

#### SHADE 99:

1lb Jar - 3738801S  
4x1lb Jar - 3738801  
5lb Jar - 3738805S  
4x5lb Jar - 3738805  
20lb Pail - 3738820

HEAT  
CURE

### REPAIR - Acrylic Denture Base

Repair Acrylic is color matched to all Garreco denture base materials. Contains special additives to prevent "yellowing" and may be used for any acrylic repair.

**INDICATION FOR USE:**  
Denture base repair

**CURE TYPE:**  
Self cure  
**INJECTABLE:** Yes  
**MICROWAVABLE:** Yes  
**LIQUID TYPE:**  
Monomer #3 Self Cure  
**POWDER/LIQUID RATIO:** 3:1



#### PINK:

1lb Jar - 3018001S  
4x1lb Jar - 3018001  
5lb Jar - 3018005S  
4x5lb Jar - 3018005  
20lb Pail - 3018020

#### PURPLE FIBER:

1lb Jar - 3013401S  
4x1lb Jar - 3013401  
5lb Jar - 3013405S  
4x5lb Jar - 3013405  
20lb Pail - 3013420

#### CLEAR:

1lb Jar - 3010101S  
4x1lb Jar - 3010101  
5lb Jar - 3010105S  
4x5lb Jar - 3010105  
20lb Pail - 3010120

#### LIGHT FIBER:

1lb Jar - 3010201S  
4x1lb Jar - 3010201  
5lb Jar - 3010205S  
4x5lb Jar - 3010205  
20lb Pail - 3010220

#### LIGHT CHARACTERIZED:

1lb Jar - 3010401S  
4x1lb Jar - 3010401  
5lb Jar - 3010405S  
4x5lb Jar - 3010405  
20lb Pail - 3010420

#### DARK FIBER:

1lb Jar - 3010301S  
4x1lb Jar - 3010301  
5lb Jar - 3010305S  
4x5lb Jar - 3010305  
20lb Pail - 3010320

#### SHADE 99:

1lb Jar - 3018801S  
4x1lb Jar - 3018801  
5lb Jar - 3018805S  
4x5lb Jar - 3018805  
20lb Pail - 3018820

SELF  
CURE

### POUR TYPE - Acrylic Denture Base

Recently reformulated, our Pour Type Acrylic has been specifically designed for use with all pour techniques. Ample pouring time and superior flexural strength are added benefits of this product.

**INDICATION FOR USE:**  
Denture base

**CURE TYPE:** Self cure  
**INJECTABLE:** No  
**MICROWAVABLE:** No  
**LIQUID TYPE:** Monomer #5  
Pour Type  
**POWDER/LIQUID RATIO:** 2:1

#### LIGHT FIBER:

1lb Jar - 3750201S  
4x1lb Jar - 3750201  
5lb Jar - 3750205S  
4x5lb Jar - 3750205  
20lb Pail - 3750220



SELF  
CURE

## AXCENT ACRYLIC™

### ORTHODONTIC - Acrylic Powder

Axcent Orthodontic Acrylic has been formulated with added stackability for orthodontic applications. Nine stock shades may be mixed to make hundreds of unique shades. Micro-fine metallic glitter may be added to any acrylic appliance for added characterization.

**INDICATION FOR USE:** Orthodontistry procedures

**CURE TYPE:** Self cure

**INJECTABLE:** No

**MICROWAVABLE:** No

**LIQUID TYPE:** Monomer #4 Hi-Clarity

**POWDER/LIQUID RATIO:** 2:1

#### CLEAR:

1lb Jar - 3760101S  
4x1lb Jar - 3760101  
5lb Jar - 3760105S  
4x5lb Jar - 3760105  
20lb Pail - 3760120



SELF  
CURE

### TRAY MATERIAL - Acrylic Designed for Custom Trays

Tray Material by Garreco is auto curing and designed for custom trays.

**INDICATION FOR USE:** Custom trays

**CURE TYPE:** Self cure

**INJECTABLE:** No

**MICROWAVABLE:** No

**LIQUID TYPE:** Monomer #6 Tray Liquid

**POWDER/LIQUID RATIO:** 2:1

#### BLUE:

1lb Jar - 3301501S  
4x1lb Jar - 3301501  
5lb Jar - 3301505S  
4x5lb Jar - 3301505  
20lb Pail - 3301520

#### WHITE:

1lb Jar - 3303001S  
4x1lb Jar - 3303001  
5lb Jar - 3303005S  
4x5lb Jar - 3303005  
20lb Pail - 3303020



SELF  
CURE

### LIQUID MONOMER

Garreco Liquid Monomers have been specifically matched to optimally perform with all Garreco powders. All liquid formulas are cross linked, "non-yellowing", and packaged in ultraviolet proof containers to ensure the longest possible shelf life.

**INDICATION FOR USE:** Liquid for acrylic denture base powders

**POWDER/LIQUID RATIO:** 2:1

#### #1 HEAT CURE:

8oz Bottle - 3780008S  
4x8oz Bottle - 3780008  
1qt Bottle - 3780030S  
4x1qt Bottle - 3780030

#### #2 15 MINUTE:

8oz Bottle - 3800008S  
4x8oz Bottle - 3800008  
1qt Bottle - 3800030S  
4x1qt Bottle - 3800030

#### #3 SELF CURE:

8oz Bottle - 3830008S  
4x8oz Bottle - 3830008  
1qt Bottle - 3830030S  
4x1qt Bottle - 3830030

#### #4 HI-CLARITY:

8oz Bottle - 3860008S  
4x8oz Bottle - 3860008  
1qt Bottle - 3860030S  
4x1qt Bottle - 3860030

#### #5 POUR TYPE:

8oz Bottle - 3910008  
4x8oz Bottle - 3910008S  
1qt Bottle - 3910030  
4x1qt Bottle - 3910030S

#### #6 TRAY LIQUID:

8oz Bottle - 3920008  
4x8oz Bottle - 3920008S  
1qt Bottle - 3920030  
4x1qt Bottle - 3920030S



#### WHICH MONOMER DO I USE?

Use the chart below to determine which monomer to use with each acrylic powder.

#### ACRYLIC POWDER:

DYNASTY HI-IMPACT  
PREMIUM  
STANDARD  
15 MINUTE  
REPAIR  
ORTHODONTIC (CLEAR)  
POUR TYPE  
TRAY MATERIAL

#### MONOMER:

#1 HEAT CURE  
#1 HEAT CURE  
#1 HEAT CURE  
#2 15 MINUTE  
#3 SELF CURE  
#4 HI-CLARITY  
#5 POUR TYPE  
#6 TRAY LIQUID

# CROWN & BRIDGE

## INVESTMENT

### ULTRAVEST™ - Investment Material

Ultravest excels when utilized with any crown & bridge alloy (except titanium) or any ceramic system. This phosphate bonded, graphite-free formula has a micro-fine particle matrix which allows precise attention to detail. This unique product may be used with either conventional or fast casting techniques and is designed for use with ring or ringless systems. UltraVest may be the only investment you will ever need!

**INDICATION FOR USE:** Investment for all crown & bridge alloys

**ITEM NUMBER:**

100x90g Pouches- 6109023



 **WARNING:** Cancer - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)

### PERFECT PRESS™ - Investment Material


Perfect Press is a graphite-free, phosphate bonded precision investment developed specially for pressed ceramics. Utilizing a fast casting technique, Perfect Press not only divests easily, but delivers accurate consistent fits. Conventional or fast casting. For use with Hi-X or Regular investment/expansion liquids.

**INDICATION FOR USE:** Investment for all crown & bridge alloys

**ITEM NUMBER:**

100x100g Pouches- 6150026



 **WARNING:** Cancer - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)

### INVESTMENT LIQUID - For Ultravest™, Perfect Press™, and Vitex™ 3D Vest

Investment Liquid for use with Ultravest, Perfect Press, and Vitex™ 3D Vest. Available in 1 liter bottles in various quantities.

**INDICATION FOR USE:** Investment material liquid

**ITEM NUMBER:**

1 Liter Bottle - 6200003S

4x1 Liter Bottle - 6200412

16x1 Liter Bottle - 6201647

ULTRAVEST LIQUID REQUIREMENTS:		PERFECT PRESS LIQUID REQUIREMENTS:		3D VEST LIQUID REQUIREMENTS:	
Case Size	Liquid	Case Size	Liquid	Case Size	Liquid
100x90g	2 Liters	100x100g	2 Liters	50x400g	3 Liters



## C&B SUPPLIES

### CRUCIBLES

Type A, non-contaminating, quartz crucibles. Available in slotted or hooded. Sold in pairs.

**INDICATION FOR USE:** Crucibles

**ITEM NUMBER:**

Hooded 2 Pack - CRU-H

Slotted 2 Pack - CRU-S





## C&B SUPPLIES

### DEBUBBLIZER - Surfactant

A superior alcohol-free surfactant for reducing surface tension on wax patterns, impressions, and models.

**INDICATION FOR USE:** Reducing surface tension on wax patterns, impressions, and models

**DEBUBBLIZER:**

8oz Bottle Blue - 3551508S  
 4x8oz Bottle Blue - 3551508  
 1qt Bottle Blue - 3551530S  
 4x1qt Bottle Blue - 3551530



### DIE SPACER - Precision Fit Agent

Die Spacer is used to provide sufficient space between the die and crown to accommodate cements used during seating. Garreco fast drying Die Spacers are formulated with additional ingredients to help seal and protect dies from wear and abrasion. Garreco's quality guarantees uniform thickness time after time. Available in 2 oz (brush cap) bottles or 4 oz (refill) bottles.

**INDICATION FOR USE:** Provides sufficient space between die and crown to accommodate cement

**BLUE SPACER:**

2oz Brush Cap Bottle - 7002002  
 12x2oz Brush Cap Bottle - 7002012  
 4oz Refill Bottle - 7002004  
 12x4oz Refill Bottle - 7002014

**GREY SPACER:**

2oz Brush Cap Bottle - 7002102  
 12x2oz Brush Cap Bottle - 7002112  
 4oz Refill Bottle - 7002104  
 12x4oz Refill Bottle - 7002114

**THINNER:**

2oz Brush Cap Bottle - 7002202  
 12x2oz Brush Cap Bottle - 7002212  
 4oz Refill Bottle - 7002204  
 12x4oz Refill Bottle - 7002214



 **WARNING:** Cancer - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)



**HAZMAT ITEM:**  
 Special shipping conditions may apply to these items.

### DISPOSABLE PLUNGERS

Universal disposable plungers for pressing ceramics. Pre-heat not required. Available in 2 sizes (2 gram and 3 gram).

**INDICATION FOR USE:** Pressing ceramics

**ITEM NUMBER:**

2g 50 Count - 6080000  
 3g 50 Count - 6080300



### FLEX TRAY

This non-stick firing tray combines the best features of rigid trays and soft liners. You'll love the flexibility, durability, and multiple uses for this tray

**INDICATION FOR USE:** Flexible tray material for firing crowns and bridges

**ITEM NUMBER:**

1 Bag (4 Trays/Bag) - 6840001



### FLEXI-RING™ - Ringless Crown & Bridge System

Advanced Ringless System allows for complete expansion of investment materials due to flexible silicone construction. Each ring is constructed with inner ribs to provide additional surface area for optimum thermal conductivity. All rings are extremely durable and designed for years of continued use in all conditions. Available in 60g, 90g, 180g, and 300g sizes.

**INDICATION FOR USE:** Ringless investment procedures

**EX-SMALL (60G):**

Ring (2.5cm/1") - 6240001  
 Base - 6260001

**SMALL (90G):**

Ring (4.3CM/1 3/4") - 6240101  
 Base - 6260101

**MEDIUM (180G):**

Ring (6CM/2 3/8") - 6240201  
 Base - 6260201

**LARGE (300G):**

Ring (7.5CM/3") - 6240301  
 Base - 6260301



## C&B SUPPLIES

### HONEYCOMB TRAYS

Specifically designed high heat firing trays for use with all furnaces and techniques. Available with ceramic or metal pins.

**INDICATION FOR USE:** Platform for firing crowns and bridges

**CERAMIC PINS:**

2 Trays with Pins - 7760801

**METAL PINS:**

2 Trays with Pins - 7760802



### VEST-IT™ - Refractory Material

Vest-it is an asbestos-free, high heat refractory material used to create custom firing pegs.

**INDICATION FOR USE:** Creating custom firing pegs

**ITEM NUMBER:**

1 Syringe (10cc) - 6840002

**VOLUME DISCOUNT ON:**

100 - 499 SYRINGES

500+ SYRINGES

Ask your dealer!



## LAB SUPPLIES

### ULTRASONIC CLEANERS

#### ULTRASONIC CLEANER #1 - General Purpose Concentrate

Ultra-Sonic Cleaner #1 is a highly concentrated formula designed to remove buffing compound, food particles, pumice, rouge, tripoli and all types of soils from dentures, bridges, crowns, orthodontic appliances, aspirator tips, burs, diamonds, glass slabs, Instruments, and more. This solution can also be used as an ultra-sonic conductor when used with other ultra-sonic cleaners in a glass beaker.

**USC #1 BLUE:**

1gal Bottle - 3570040S

4x1gal Bottle - 3570040



#### ULTRASONIC CLEANER #2 - Tartar and Stain Remover

Ultra-Sonic Cleaner #2 is specially formulated to easily remove calculus, food stains, permanent cement, tartar and tobacco from dentures, bridges, crowns, orthodontic appliances, and other dental restorations.

**USC #2 CLEAR:**

1gal Bottle - 3590040S

4x1gal Bottle - 3590040



#### ULTRASONIC CLEANER #5 - Plaster and Stone Remover

Ultra-Sonic Cleaner #5 is especially formulated to remove plaster and stone from dentures, bridges, crowns, orthodontic appliances, and other dental restorations.

**USC #5 PINK:**

1gal Bottle - 3610040S

4x1gal Bottle - 3610040



## GENERAL USE

### SEPARATING FILM

Garreco Separating Film is a superior alternative to tinfoil and rinses completely from teeth and metal frameworks to provide unparalleled separation from the acrylic and mold. This product is available in Regular and Concentrate.

**INDICATION FOR USE:** Separating film

#### REGULAR RED:

1qt Bottle - 3458330S  
4x1qt Bottle - 3458330  
1gal Bottle - 3458340S  
4x1gal Bottle - 3458340

#### REGULAR CLEAR:

1qt Bottle - 3450130S  
4x1qt Bottle - 3450130  
1gal Bottle - 3450140S  
4x1gal Bottle - 3450140

#### CONCENTRATE RED:

1gal Bottle - 3468330S  
4x1gal Bottle - 3468330



**WARNING:** Cancer - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)

## PERFORM™ PRODUCTS


### PERFORM™ TOOTH SIEVE

PERform Tooth Sieves are a must when using pour techniques. These compartmented plastic sieves are heat resistant and durable for years of service.

**INDICATION FOR USE:** Used for easy scalding of teeth

**ITEM NUMBER:**

PERform Tooth Sieve Blue - 4100130

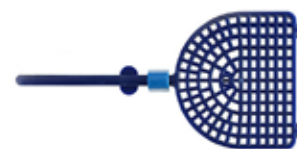

### PERFORM™ TOOTH CUTTING BUR

PERform Tooth Cutting Burs are used for cutting diatorics in acrylic teeth.

**INDICATION FOR USE:** Cutting diatorics

**ITEM NUMBER:**

PERform Tooth Cutting Bur - 4100040



### PERFORM™ FILTER PAPERS

PERform™ Filter Papers for use with PERform™ unit.

**INDICATION FOR USE:** Filter papers

**ITEM NUMBER:**

PERform Filter Papers 500/pack - 4100100



# PARTIAL DENTURES


## VITEX™ PRODUCTS


### VITEX™ INVESTMENT


Vitex (Ethyl Silicate) partial denture framework investment material has been developed especially for producing chrome cobalt castings. This product offers precise details due to an ultra fine particle size and unique two part blend of liquids. May be used for all high heat techniques. Formulated to mix and pour easily, Vitex investment maintains excellent integrity when used in conjunction with Vitex Model Hardener.

**INDICATION FOR USE:** Investment material

#### 100LB KIT:

100lb Kit (3 Cartons) - 5007000   
Includes: 2 x 50lb Investment Sand, 2 x 118oz "A" Liquid, and 2 x 14.8oz "B" Liquid

100lb Kit Liquid Only (Carton) - 5012200   
Includes: 2 x 118oz "A" Liquid, and 2 x 14.8oz "B" Liquid

50lb Sand Only (Carton) - 5008000   
118oz "A" Liquid (Carton) - 5009000  
14.8oz "B" Liquid (Carton) - 5011000

#### 50LB KIT:

50lb Kit (2 Cartons) - 5006000  
Includes: 1 x 50lb Investment Sand, 4 x 28oz "A" Liquid, and 4 x 3.7oz "B" Liquid

28oz "A" Liquid (Carton) - 5010000  
3.7oz "B" Liquid (Carton) - 5012000



 Denotes an item that cannot be shipped by FedEx Ground; must be shipped FedEx Freight only



**HAZMAT ITEM:**  
Special shipping conditions may apply to these items.

**WARNING:** Cancer - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)

## VITEX™ PRODUCTS

### VITEX™ 3D VEST

Vitex™ 3D Vest is a phosphate-bonded, graphite free and very fine grain precision investment for partial denture frameworks made with “Rapid-Prototyping” technologies and is specially developed to burn out plastic frames made with Rapid-Prototyping resins. 3D Vest is usable for both the conventional and speed casting techniques and is suitable for casting all types of dental alloys (except Titanium).

**INDICATION FOR USE:** Investment material designed for 3D Printed/Milled Wax and Plastic Scaffolding.

**ITEM NUMBER:**

50x400g Pouches - 5013400

USE WITH INVESTMENT LIQUID ON PAGE 16

 **WARNING:** Cancer - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)



### VITEX™ REFRACTORY MODEL HARDENER

Vitex Refractory Model Hardener (“Cold Dip”) is designed to be used at room temperature with a warm model. Use of this material exhibits the same results as older techniques, while saving time and money.


**INDICATION FOR USE:** Model hardener

**ITEM NUMBER:**

1qt Bottle - 5005030S

4x1qt Bottle - 5005030

1gal Bottle - 5005000 

4x1gal Bottle - 5005040 



**HAZMAT ITEM:** Special shipping conditions may apply to this item.



Denotes an item that cannot be shipped by FedEx Ground; must be shipped FedEx Freight only

### VITEX™ DUPLICATING MATERIAL

Although developed for use in the construction of partial denture frameworks, this versatile hydrocolloid material has proven to be excellent for a variety of dental related tasks. With superior reproduction of detail and endurance for multiple melting cycles. Available in two types (all purpose and regular). All Purpose melts at 190° and holds at 135°. Regular melts at 200° and holds at 125°. Available in Neutral or Blue.

**INDICATION FOR USE:** Partial denture framework construction

**ALL PURPOSE (Melt at 190°F, hold at 135°F):**

2gal Blue - 5001015

5gal Blue - 5001050

2gal Neutral - 5000015

5gal Neutral - 5000050

**REGULAR (Melt at 200°F, hold at 125°F)**

2gal Blue - 5001515

5gal Blue - 5001550



### VITEX™ ALLOY

This unique blend of metals provides short term flexibility and long term rigidity required by today's advanced technicians. Vitex partial denture alloy has been developed for use with all high heat techniques, and is priced well below other premium alloys.

**INDICATION FOR USE:** Dental alloy

**ITEM NUMBER:**

1lb Alloy Bullets Carton - 5018500



## VITEX™ PRODUCTS

### VITEX™ SUPERFORMS

Vitex Superform Patterns are made from a wax/plastic blend. This mix offers all the advantages of both components with optimum burnout. Available in multiple styles in protective packaging. 10 cards per box.

**INDICATION FOR USE:** Patterns for casting alloy



**B-1**  
10 CARD BOX - 5051010



**B-2**  
10 CARD BOX - 5052010



**B-4**  
10 CARD BOX - 5054010



**PI-PE**  
10 CARD BOX - 5113010



**G-45/G-55**  
10 CARD BOX - 5083010  
10 CARD BOX - 5084010



**GC-1**  
10 CARD BOX - 5085010



**RC-1**  
10 CARD BOX - 5103510



**GR-4**  
10 CARD BOX - 5099010



**RI-2**  
10 CARD BOX - 5105010



**RR-3**  
10 CARD BOX - 5108010



**GC-3**  
10 CARD BOX - 5087010



**GD-1**  
10 CARD BOX - 5093010



**GR-5**  
10 CARD BOX - 5101010



**BR-1**  
10 CARD BOX - 5072010



**BP-1**  
10 CARD BOX - 5068010



**GC-5**  
10 CARD BOX - 5089010



**GD-2**  
10 CARD BOX - 5094010



**GR-6**  
10 CARD BOX - 5102010



**BC-1**  
10 CARD BOX - 5058010



**BP-11**  
10 CARD BOX - 5071010



**GC-6**  
10 CARD BOX - 5091010



**GD-3**  
10 CARD BOX - 5095010



**GR-6**  
10 CARD BOX - 5102010



**RT-0**  
10 CARD BOX - 5109010



**BL-1**  
10 CARD BOX - 5063010



**GC-7**  
10 CARD BOX - 5092010




**GD-4**  
10 CARD BOX - 5097010



**CR-1**  
10 CARD BOX - 5076010



**RT-1**  
10 CARD BOX - 5111010

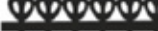

**SI-1**  
10 CARD BOX - 5112010



**GD-3L**  
10 CARD BOX - 5096010



**GD-5**  
10 CARD BOX - 5098010



**RI-1**  
10 CARD BOX - 5104010



**FC-0**  
10 CARD BOX - 5079010

### SUBSCRIBE AND SAVE

Join our mailing list to get special sales and promotions delivered directly to your inbox.  
Visit [Garreco.com/signup](http://Garreco.com/signup) today to get started. Check out [Garreco.com/promotions](http://Garreco.com/promotions) to see the current promotions we have available.



[Garreco.com/signup](http://Garreco.com/signup)



[Garreco.com/promotions](http://Garreco.com/promotions)



# ROTARY

## PERLADIA POLISHERS

### PERLADIA POLISHERS - Diamond Impregnated and Long Lasting



Perlada Polishers have been the standard in porcelain polishing for many years and can be used with any porcelain. They are long lasting, generate no heat, and eliminate the need for re-glazing. No porcelain chipping and no polishing paste. All Perlada Polishers are impregnated with diamonds!

**INDICATION FOR USE:** Shaping, finishing and polishing with no chipping or heat generation.

#### 25MM - KNIFE EDGE:




Ochre  
Extra Coarse Gross  
Trim Zirconia  
PDS25XCZO



Red  
Coarse  
Quick Reduction  
PDS25CR



Tan  
Medium- Coarse  
Contouring  
PDS25MCT

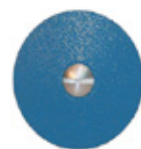

Blue  
Medium  
Finishing  
PDS25MB



Green  
Fine  
Initial Polish  
PDS25FG



Purple  
Extra Fine  
Pre-Polish  
PDS25XFP



Yellow  
High Shine  
Glaze  
PDS25HSY



Black  
Coarse  
Zirconia  
PDS25CZBK

## NATURALGLAZE DIAMOND POLISHERS

### NATURALGLAZE - Diamond Polishing System

Garreco ensures lasting good looks and the highest quality for all your restorations with Garreco rotary instruments. Naturalglaze polishers are excellent for shaping, finishing, and polishing ceramic materials with no chipping or heat generation.

**INDICATION FOR USE:** Shaping, finishing, and polishing ceramic with no chipping or heat generation.

#### 25MM - KNIFE EDGE:




Coarse  
Solid Disc  
U6425-250-1



Medium  
Solid Disc  
U6325-250-1



Fine  
Solid Disc  
U6225-250-1

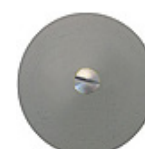

High Shine  
Solid Disc  
U6125-250-1



Coarse  
Plastic Center  
6425-250-1



Medium  
Plastic Center  
6325-250-1



Fine  
Plastic Center  
6225-250-1



High Shine  
Plastic Center  
6125-250-1

#### 17MM - SQUARE EDGE:



Coarse  
Plastic Center  
6409-170-1



Medium  
Plastic Center  
6309-170-1



Fine  
Plastic Center  
6209-170-1



High Shine  
Plastic Center  
6109-170-1

#### FLAME POLISHERS:



Medium  
6360-1



Fine  
6260-1



High Shine  
6160-1

## CARBIDES - FG & HP

### FG CARBIDES (Friction Grip)

Garreco ensures lasting good looks and the highest quality for all your restorations with Garreco rotary instruments.

**INDICATION FOR USE:** Shaping, finishing, and polishing ceramic with no chipping or heat generation.

#### FG CARBIDES:

FG 1 CARBIDE x 10 - FG1	FG 33 1/2 CARBIDE x 10 - FG33 1/2	FG 558 CARBIDE x 10 - FG558
FG 1/2 CARBIDE x 10 - FG1/2	FG 330 CARBIDE x 10 - FG330	FG 57 CARBIDE x 10 - FG57
FG 1/4 CARBIDE x 10 - FG1/4	FG 331 CARBIDE x 10 - FG331	FG 58 CARBIDE x 10 - FG58
FG 1156 CARBIDE x 10 - FG1156	FG 34 CARBIDE x 10 - FG34	FG 6 CARBIDE x 10 - FG6
FG 1157 CARBIDE x 10 - FG1157	FG 35 CARBIDE x 10 - FG35	FG 700 CARBIDE x 10 - FG700
FG 1158 CARBIDE x 10 - FG1158	FG 36 CARBIDE x 10 - FG36	FG 701 CARBIDE x 10 - FG701
FG 1557 CARBIDE x 10 - FG1557	FG 37 CARBIDE x 10 - FG37	FG 701L CARBIDE x 10 - FG701L
FG 1558 CARBIDE x 10 - FG1558	FG 39 CARBIDE x 10 - FG39	FG 702 CARBIDE x 10 - FG702
FG 1558 CARBIDE x 100 - FG1558X100	FG 4 CARBIDE x 10 - FG4	FG 703 CARBIDE x 10 - FG703
FG 169L CARBIDE x 10 - FG169L	FG 556 CARBIDE x 10 - FG556	7901 FINISHING CARBIDE 12 BLADES x 10 - FG7901
FG 171 CARBIDE x 10 - FG171	FG 556L CARBIDE x 10 - FG556L	7903 FINISHING CARBIDE 12 BLADES x 10 - FG7903
FG 2 CARBIDE x 10 - FG2	FG 557 CARBIDE x 10 - FG557	FG 8 CARBIDE x 10 - FG8
FG 331 CARBIDE x 10 - FG331	FG 557L CARBIDE x 10 - FG557L	FG METAL/CROWN CUTTING CARBIDE x10 - FGPR2



### HP CARBIDES (Hand Piece)

Garreco ensures lasting good looks and the highest quality for all your restorations with Garreco rotary instruments.

**INDICATION FOR USE:** Shaping, finishing, and polishing ceramic with no chipping or heat generation.

#### HP CARBIDES:

HP 1 STRAIGHT CARBIDE x 10 - HP1	HP 35 STRAIGHT CARBIDE x 10 - HP35	HP 703 STRAIGHT CARBIDE x 10 - HP703
HP 1/2 STRAIGHT CARBIDE x 10 - HP1/2	HP 36 STRAIGHT CARBIDE x 10 - HP36	HP 8 STRAIGHT CARBIDE x 10 - HP8
HP 1/4 STRAIGHT CARBIDE x 10 - HP1/4	HP 37 STRAIGHT CARBIDE x 10 - HP37	RA6 RIGHT ANGLE LATCH
HP 10 STRAIGHT CARBIDE x 10 - HP10	HP 4 STRAIGHT CARBIDE x 10 - HP4	CARBIDE x 10 - RA6
HP 1157 STRAIGHT CARBIDE x 10 - HP1157	HP 557 STRAIGHT CARBIDE x 10 - HP557	2 RIGHT ANGLE LATCH SURGICAL
HP 1158 STRAIGHT CARBIDE x 10 - HP1158	HP 558 STRAIGHT CARBIDE x 10 - HP558	CARBIDE x 10 - RAS2
STRAIGHT CARBIDE x 10 - HP 1556	HP 57 STRAIGHT CARBIDE x 10 - HP57	
HP 1557 STRAIGHT CARBIDE x 10 - HP1557	HP 6 STRAIGHT CARBIDE x 10 - HP6	
HP 1558 STRAIGHT CARBIDE x 10 - HP1558	HP 699 STRAIGHT CARBIDE x 10 - HP699	
HP 2 STRAIGHT CARBIDE x 10 - HP2	HP 7 STRAIGHT CARBIDE x 10 - HP7	
HP 3 STRAIGHT CARBIDE x 10 - HP3	HP 701 STRAIGHT CARBIDE x 10 - HP701	
HP 34 STRAIGHT CARBIDE x 10 - HP34	HP 701L STRAIGHT CARBIDE x 10 - HP701L	



## HUB BRUSHES

### PLASTIC & WOOD HUB BRUSHES

Garreco ensures lasting good looks and the highest quality for all your restorations with Garreco rotary instruments.

**INDICATION FOR USE:** Hub brushes for shaping, finishing, and polishing.

#### WOOD HUB BRUSHES:

WOOD HUB BRUSH WHEEL B-12 (12 PACK) - 734B122  
 WOOD HUB BRUSH WHEEL B-20 (12 PACK) - 734B242  
 WOOD HUB BRUSH WHEEL B-24 (12 PACK) - 734B202  
 WOOD HUB BRUSH WHEEL B-25 (12 PACK) - 734B252  
 WOOD HUB BRUSH WHEEL B-27 (12 PACK) - 734B272  
 WOOD HUB BRUSH WHEEL B-29 (12 PACK) - 734B292

#### PLASTIC HUB BRUSHES:

PLASTIC HUB BRUSHES B-12 (12 PACK) - 735B122  
 PLASTIC HUB BRUSHES B-14 (12 PACK) - 735B142  
 PLASTIC HUB BRUSHES B-20 (12 PACK) - 735B202  
 PLASTIC HUB BRUSHES B-24 (12 PACK) - 735B242



## MUSLIN BUFFS & BRISTLE BRUSHES


### MUSLIN BUFF WHEELS

Garreco high-quality muslin buff wheels for polishing and finishing available in several varieties.

**INDICATION FOR USE:** Muslin buffs for polishing and finishing.

#### MUSLIN BUFF WHEEL - SHELLAC CENTER:

MUSLIN BUFF WHEEL SHELLAC CENTER (4X30) (12 PAK) - 7340130  
 MUSLIN BUFF WHEEL SHELLAC CENTER (4X42) (12 PAK) - 7340142  
 MUSLIN BUFF WHEEL SHELLAC CENTER (2 X16) (12 PAK) - 7340216  
 MUSLIN BUFF WHEEL SHELLAC CENTER (2X30) (12 PAK) - 7340230  
 MUSLIN BUFF WHEEL SHELLAC CENTER (3X30) (12 PAK) - 7340330  
 MUSLIN BUFF WHEEL SHELLAC CENTER (3X35) (12 PAK) - 7340335  
 MUSLIN BUFF WHEEL SHELLAC CENTER (4X24) (12 PAK) - 7340424  
 MUSLIN BUFF WHEEL SHELLAC CENTER (5X30) (12 PAK) - 7340530  
 MUSLIN BUFF WHEEL SHELLAC CENTER (5X45) (12 PAK) - 7340545  
 MUSLIN BUFF WHEEL SHELLAC CENTER (6X30) (12 PAK) - 7340630

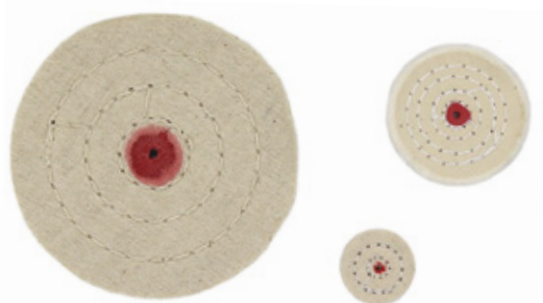


#### MUSLIN BUFF WHEEL:

MUSLIN BUFF WHEEL (1X16) (12 PAK) - 7340116

#### FINEX MUSLIN BUFF WHEEL - PLASTIC CENTER:

FINEX MUSLIN BUFF WHEEL UN-STITCHED PLASTIC CENTER (4X30) (12 PAK) - 7344412

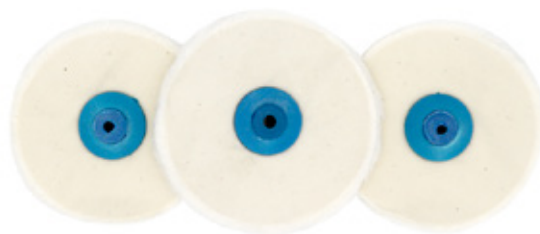

### BRISTLE BRUSHES

Garreco ensures lasting good looks and the highest quality for all your restorations with Garreco rotary instruments.

**INDICATION FOR USE:** Shaping, finishing, and polishing.

#### BRISTLE BRUSHES #9:


BRISTLE BRUSHES #9 SOFT (9/16") (12 PAK) - 7500001  
 BRISTLE BRUSHES #9 SOFT (9/16") (12X12 PAK) - 7500002  
 BRISTLE BRUSHES #9 STIFF (9/16") (12 PAK) - 7520001  
 BRISTLE BRUSHES #9 STIFF (9/16") (12X12 PAK) - 7520002  
 BRISTLE BRUSHES #9 X-STIFF (9/16") (12 PAK) - 7540001  
 BRISTLE BRUSHES #9 X-STIFF (9/16") (12X12 PAK) - 7540002

#### BRISTLE BRUSHES #12:

BRISTLE BRUSHES #12 SOFT (3/4") (12 PAK) - 7500005  
 BRISTLE BRUSHES #12 SOFT (3/4") (12X12 PAK) - 7500006  
 BRISTLE BRUSHES #12 STIFF (3/4") (12 PAK) - 7520005  
 BRISTLE BRUSHES #12 STIFF (3/4") (12X12 PAK) - 7520006  
 BRISTLE BRUSHES #12 X-STIFF (3/4") (12 PAK) - 7540005  
 BRISTLE BRUSHES #12 X-STIFF (3/4") (12X12 PAK) - 7540006

#### BRISTLE BRUSHES #11:

BRISTLE BRUSHES #11 SOFT (11/16") (12 PAK) - 7500003  
 BRISTLE BRUSHES #11 SOFT (11/16") (12X12 PAK) - 7500004  
 BRISTLE BRUSHES #11 STIFF (11/16") (12 PAK) - 7520003  
 BRISTLE BRUSHES #11 STIFF (11/16") (12X12 PAK) - 7520004  
 BRISTLE BRUSHES #11 X-STIFF (11/16") (12 PAK) - 7540003  
 BRISTLE BRUSHES #11 X-STIFF (11/16") (12X12 PAK) - 7540004

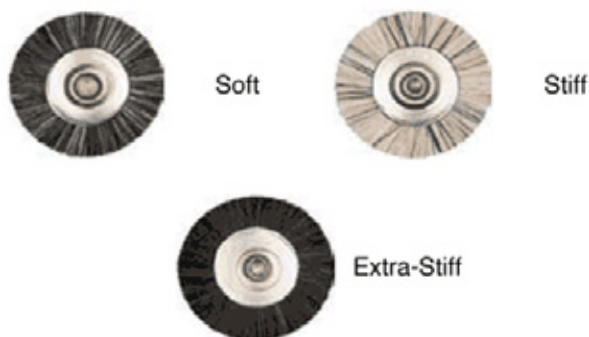

## INTERPROXIMAL ROTARY

Garreco has several dental rotary discs that are specifically designed to reach interproximal areas for Ceramics and Zirconia/Pressables.

**INDICATION FOR USE:** Shaping, finishing and polishing interproximal areas.

**20MM - INTERPROXIMAL:**



NATURALGLAZE  
Medium  
Ceramics  
SL20DG




NATURALGLAZE  
Fine  
Ceramics  
SL20DMF



NATURALGLAZE  
High Shine  
Ceramics  
SL20DG



PERLADIA S-20-M  
Medium  
Zirconia/Pressables  
6125-200-1M

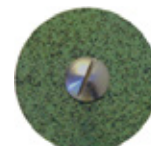

POLIDIA ULTRA-THIN  
Medium  
Zirconia/Pressables  
6113-200-1

## SILICONES

### SILICONES

#### SILAFORM™ - Silicone Putty

Silaform silicone putty is a one to one mix material which is excellent for a variety of procedures. Silaform 90 (blue) has a Shore A hardness of 90 while Silaform 70 (purple) is more flexible and has a Shore A hardness of 70. Both materials resist deformation and abrasion normally associated with silicone putty. The absence of expansion and contraction offers endless possibilities for this economically priced product from Garreco. Silaform is only available in sets.

**INDICATION FOR USE:** Duplicating material

##### SILAFORM™ 70 HARD (PURPLE):

1kg Set (Canister) - 6803404  
25kg Set (Pail) - 6803455

##### SILAFORM™ 90 EXTRA HARD (BLUE):

1.5kg Set (Canister) - 6801505  
25kg Set (Pail) - 6801555



#### ACCUSIL™ - Silicone Duplicating Material

A high quality platinum based silicone duplicating material for use with all techniques. Accuracy is constant and predictable at 0.0% expansion/contraction. Excellent strength of 16-18 Shore A hardness. May be poured multiple times with gypsum or investment without deterioration. Suitable for use in dispensers. Accusil is only available in sets.

**INDICATION FOR USE:**  
Duplicating material

##### ACCUSIL:

1kg Set (1 Activator & 1 Base) - 6705005  
25kg Set (1 Activator & 1 Base) - 6705055

**MIXING RATIO:** 1:1, A/B  
**WORKING TIME:** >5 min  
**CURING TIME:** 30-45 min  
**ELONGATION OF TEAR:** ~310%



## SILICONES

### RAPIDSIL™ - Silicone Duplicating Material

Lower cost alternative to Accusil with a slightly reduced Shore A hardness of 15 and a faster curing time. Rapidsil offers no expansion or contraction. Due to the use of a platinum base, this material may be poured hundreds of times with equally pleasing results. Rapidsil is only available in sets.

**INDICATION FOR USE:** Duplicating material

**MIXING RATIO:** 1:1, A/B

**WORKING TIME:** >4 min

**CURING TIME:** 15-20 min

**ELONGATION OF TEAR:** ~300%

**RAPIDSIL:**

1kg Set (1 Activator & 1 Base) - 6701805

25kg Set (1 Activator & 1 Base) - 6701855



### SILSAVER™

Especially designed plastic spacers for silicone duplication. Adds stability to silicone, while saving material cost. Reusable and available in 3 sizes.

**INDICATION FOR USE:** Spacers for silicone duplication

**SILSAVERS:**

Small (each) - 6480101

Medium (each) - 6480201

Large (each) - 6480301



### DUPLICATING CROSS - Placement Device

Re-useable, threaded placement device for Silsaunders. Item discontinued; limited quantity available.

**INDICATION FOR USE:** Silsaver placement

**ITEM NUMBER:**

Duplicating Cross - 6440003



### SPRUE FORMERS

Specially designed plastic sprue formers for use with refractory casting techniques.

**INDICATION FOR USE:** Refractory casting

**ITEM NUMBER:**

Sprue Formers (10/pkg) - 6600001



### DUPLICATING RING

Disposable, clear, self adhesive duplicating ring, for use with silicone duplication techniques.

**INDICATION FOR USE:** Duplication techniques

**ITEM NUMBER:**

Duplicating Ring (40m/Roll) - 6460002





# WAXES

## BASE PLATE

### BASE PLATE WAX

Our dependable Base Plate Wax is available in Medium Soft, Medium Hard, or All Season. Medium Hard is formulated with higher melting waxes and is designed to work in warmer climates. Medium Soft is designed to be used most anywhere a good quality wax is needed. Garreco Base Plate Wax is dependable, boils out clean, is competitively priced, and polishes to a high shine.

**Sizes:** 76mm x 152mm x .060 thick, available in 1lb or 5lbs

#### MEDIUM SOFT:

1lb Box - 51510101  
5lb Box - 51510105

#### MEDIUM HARD:

1lb Box - 51510201  
5lb Box - 51510205

#### ALL SEASON:

1lb Box - 51510001  
5lb Box - 51510005



## ORTHODONTIC

### ORTHOWAX COMPACTS

Our long-lasting OrthoWax is recommended for relief from irritation that may be caused by wires and brackets and is FDA cleared. OrthoWax is available in recloseable clear plastic boxes.

OrthoWax Compacts Braces Relief Wax - 50 or 100 compacts per box.

**Available Scents:** Cinnamon, Wintergreen, Cherry, Bubblegum, and Unscented.

#### ORTHOWAX COMPACTS:

Unscented 50 Pack Box - 56720050  
Unscented - 100 Pack Box 56720099

Wintergreen 50 Pack Box - 56722750  
Wintergreen 100 Pack Box - 56722799

Bubble Gum 50 Pack Box - 56728650  
Bubble Gum 100 Pack Box - 56728699

Cherry 50 Pack Box - 56728450  
Cherry 100 Pack Box - 56728499

Cinnamon 50 Pack Box - 56728250  
Cinnamon 100 Pack Box - 56728299



## FIND A GARRECO DEALER NEAR YOU

As a manufacturer of premium dental laboratory products, we are proud of our strong dealership network that spans multiple countries across the globe. Based on your location, one of our sales experts can provide you with a complete list of authorized Garreco dealers in your area.

To find an authorized Garreco Dealer near you, please call **1-800-334-1443** or email us at [sales@garreco.com](mailto:sales@garreco.com)



## UTILITY

### REGULAR BOXING WAX

A pliable wax that does not need heating yet is easy to work with and is strong and tough. The wax is designed to reduce the possibility of distorting the peripheral roll when boxing an alginate or hydrocolloid impression.

#### BOXING WAX:

White Box - 51510201

Red Box - 51510205



### SQUARE ROPES

A soft, pliable, tacky wax that can be used for boxing alginates and hydrocolloids and building up impression trays. It is sufficiently soft to adapt itself to allow the tray to be seated. It will adhere to most any material and can be used for post damming and boxing.

#### SQUARE ROPES:

White Box - 53410000

Red Box - 53418300



### LARGE UTILITY STRIPS

Large Utility Wax Strips are sticky enough to adhere without softening, but not too sticky to be removed cleanly. These 11" long wax strips are perfect for border molding impression trays and beading impressions to preserve peripheral roll.

#### LARGE UTILITY STRIPS:

White Strips Box - 53710000

Red Strips Box - 53718300



### SMALL UTILITY STRIPS

Small Utility Wax Strips are sticky enough to adhere without softening, but not too sticky to be removed cleanly. These 11" long wax strips are perfect for border molding impression trays and beading impressions to preserve peripheral roll.

#### SMALL UTILITY STRIPS:

White Strips Box - 54710000

Red Strips Box - 54718300




### PERIPHERY STRIPS

Periphery Wax can be used with all types of impression materials for making wax stops, post-damming and lining of the periphery trays.

#### PERIPHERY STRIPS:

White Strips Box - 58719900

Blue Strips Box - 58711500

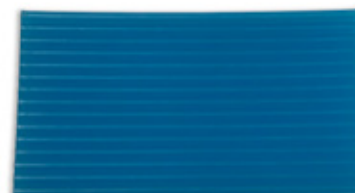

## CONTACT US



**Normal Operating Hours:**  
Monday-Thursday, 7:00am to 5:30pm



**Physical Address:**  
Garreco, LLC  
430 Hiram Road  
Heber Springs, AR 72543



**Mailing Address:**  
Garreco, LLC  
PO Box 1258  
Heber Springs, AR 72543



**Toll Free Number:** (800) 334-1443  
**Local Number:** (501) 362-6261  
**Fax Number:** (501) 362-2264



**Email Addresses:**  
**General Info:** [info@garreco.com](mailto:info@garreco.com)  
**Sales:** [sales@garreco.com](mailto:sales@garreco.com)  
**Support:** [support@garreco.com](mailto:support@garreco.com)  
**Orders:** [orders@garreco.com](mailto:orders@garreco.com)



**Online:**  
**Websites:**  
[Garreco.com](http://Garreco.com)  
[Garrecoprint.com](http://Garrecoprint.com)

**Social Media:**  
**Facebook** - [Facebook.com/GarrecoMfg](https://Facebook.com/GarrecoMfg)  
**Instagram** - [Instagram.com/garrecodental](https://Instagram.com/garrecodental)  
**YouTube** - [Youtube.com/channel/UCHZtuMaiGnf6OEzgC8TXICQ](https://Youtube.com/channel/UCHZtuMaiGnf6OEzgC8TXICQ)  
**Twitter/X** - [Twitter.com/GarrecoDental](https://Twitter.com/GarrecoDental)  
**LinkedIn** - [Linkedin.com/company/garreco-dental](https://Linkedin.com/company/garreco-dental)